

## Mouse GDF8 / Myostatin protein, His tag

Cat. No. GTX00308-pro

实验应用 Functional Assay

Speciess Mouse

Package

10 µg

## 实验应用

## 应用说明

Myostatin (MSTN) also known as growth differentiation factor 8 (GDF-8) a myokine, a protein produced and released by myocytes. It inhibit myogenesis including muscle cell growth and differentiation. Myostatin is a secreted growth differentiation factor that is a member of the TGF beta protein family. Besides, Follistatin Like Protein 3 (FSTL3) has been identified as an interactor of MSTN, thus a binding ELISA assay was conducted to detect the interaction of recombinant mouse MSTN and recombinant mouse FSTL3. Briefly, MSTN were diluted serially in PBS, with 0.01% BSA (pH 7.4). Duplicate samples of 100 µl were then transferred to FSTL3-coated microtiter wells and incubated for 2h at 37°C. Wells were washed with PBST and incubated for 1h with anti-MSTN pAb, then aspirated and washed 3 times. After incubation with HRP labelled secondary antibody, wells were aspirated and washed 3 times. With the addition of substrate solution, wells were incubated 15-25 minutes at 37°C. Finally, add 50 µl stop solution to the wells and read at 450nm immediately. The binding activity of MSTN and FSTL3 was in a dose dependent manner.

Observed MW 15 kDa.

## 属性

形式 Lyophilized powder

存储溶液 Reconstitute with 20mM Tris and 150mM NaCl to 0.1-1.0mg/ml. Do not vortex. Lyophilized from 20mM Tris, 150mM NaCl, 1mM EDTA, 1mM DTT, 0.01% SKL, 5% Trehalose.

保存剂 ProClin 300

存放说明 For short-term storage (1-2 weeks), store at 4°C. For long-term storage, store at -20°C or below. After reconstitution, keep as concentrated solution. Avoid freeze-thaw cycles.

Region/Sequence N-terminal His-Tag; Asp268~Ser376 (NP\_034964.1)

表达系统 E. coli

纯度 &gt; 95%

内毒素 &lt; 1 EU/µg


偶联 Unconjugated

注意事项 For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.




For full product information, images and publications, please visit our [website](#).

## 產品圖片

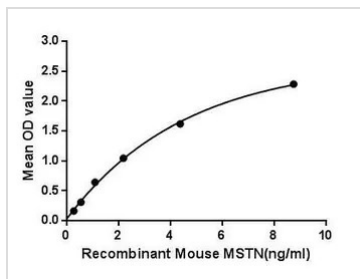


### GTx00308-pro Functional Assay Image

Functional ELISA analysis of GTx00308-pro Mouse GDF8 / Myostatin protein which can bind immobilized FSTL3 protein.

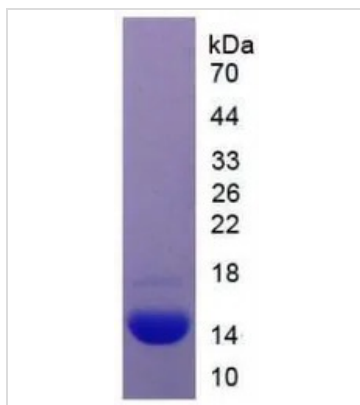

### GTx00308-pro Image

SDS-PAGE analysis of GTx00308-pro Mouse GDF8 / Myostatin protein.

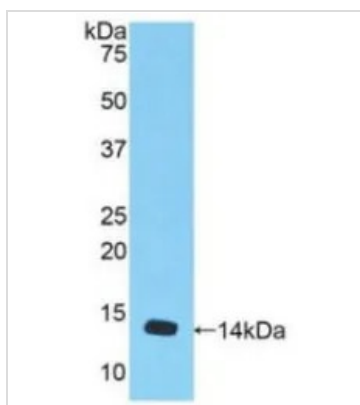

### GTx00308-pro Image

WB analysis of GTx00308-pro Mouse GDF8 / Myostatin protein.



For full product information, images and publications, please visit our [website](https://www.genetex.com).