

TER-119 antibody [TER-119] (PE)

Cat. No. GTX01475-08

宿主	Rat
克隆	Monoclonal
同种型	IgG2b
实验应用	FCM
种属反应	Mouse

引用文献 (7)

产品说明

摘要

The TER-119 antibody is named for the antigen to which it binds, a 52 kDa surface protein that is associated with glycophorin-A. TER-119 is considered to be a lineage marker for later stages of erythroid cell development, as its expression begins at the pro-erythroblast stage. TER-119 antigen is not expressed at either BFU-E or CFU-E stages, i.e. prior to the pro-erythroblast stage.

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
FCM	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹
ICC/IF: 细胞染色
IHC-P: 石蜡切片
IHC-Fr: 冰冻切片

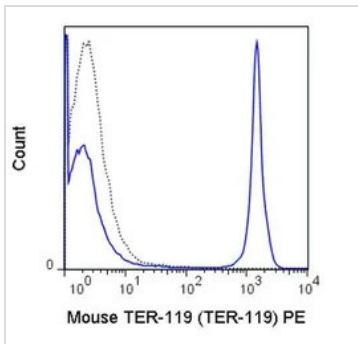
属性

形式	Liquid
存储溶液	10mM NaH ₂ PO ₄ , 150mM NaCl, 0.1% Gelatin
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. Store at 4°C. DO NOT FREEZE. Protect from light.
浓度	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Phycoerythrin (PE) Wavelength
注意事项	仅供实验室使用。不适用于人类或动物的任何临床，治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

**GTX01475-08 FCM Image**

FACS analysis of mouse C57Bl/6 bone marrow cells using GTX01475-08 TER-119 antibody [TER-119] (PE).

Solid line : primary antibody

Dashed line : isotype control

antibody amount : 0.5 μ g (5 μ l)



For full product information, images and publications, please visit our [website](#).