

Parathyroid Hormone antibody [105G7]

Cat. No. GTX01891

宿主	Mouse
克隆	Monoclonal
同种型	IgG2a
实验应用	IHC-P
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IHC-P	1:300

Note : The recommended fixative is 10% neutral-buffered formalin for paraffin-embedded tissue sections. Heat Induced Epitope Retrieval (HIER) was performed in Epitope Retrieval Solution pH6.

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	Tissue culture supernatant
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. Store at 4°C. DO NOT FREEZE.
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

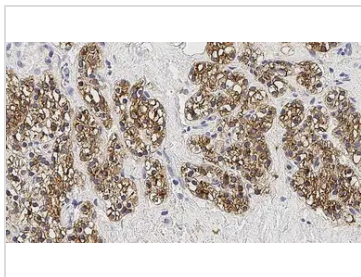

GTx01891 IHC-P Image

IHC-P analysis of human parathyroid adenoma tissue using GTx01891 Parathyroid Hormone antibody [105G7]. Note membrane/cytoplasmic staining of tumor cells.



For full product information, images and publications, please visit our [website](#).