

Collagen IV antibody [COL4/4241R]

Cat. No. GTX02748

宿主	Rabbit
克隆	Monoclonal
同种型	IgG
实验应用	IHC-P
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IHC-P	1-2 µg/ml

Note : Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris buffer with 1mM EDTA (pH 9.0) for 45 min at 95°C followed by cooling at RT for 20 minutes.

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 0.05% BSA
保存剂	0.05% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. Store at 4°C.
浓度	200 µg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

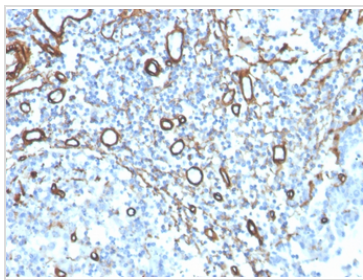

GTX02748 IHC-P Image

IHC-P analysis of human tonsil tissue section using GTX02748 Collagen IV antibody [COL4/4241R].

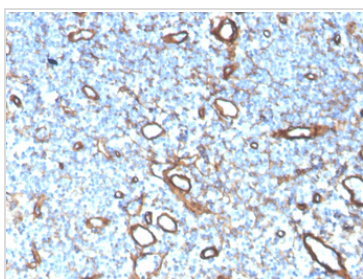

GTX02748 IHC-P Image

IHC-P analysis of human tonsil tissue section using GTX02748 Collagen IV antibody [COL4/4241R].



For full product information, images and publications, please visit our [website](#).