

IL-12 p70 ELISA pair [MT86/221/MT704]

Cat. No. GTX03071

实验应用	ELISA, Sandwich ELISA
种属反应	Human

Package

1 pair

产品说明

Product Content

Cat No.	Type	宿主	偶联	Package
GTX02957	Capture antibody	Mouse	Unconjugated	250 µg
GTX03013-02	Detection antibody	Mouse	Biotin	250 µg

实验应用

产品说明

GTX02957 (MT86/MT221) is able to detect native and recombinant human p40 subunit of IL-12 and IL-23. In combination with GTX03013-02 (MT704), it recognizes IL-12 p70 heterodimer (p40 + p35).

属性

形式	Liquid
----	--------


存放说明

Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.

注意事项

仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

产品图片

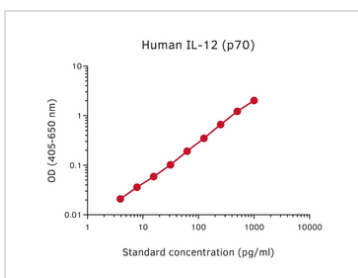

GTX03071 ELISA Image

Sandwich ELISA analysis of human IL-12 (p70) protein using GTX02957 IL12B / IL12 p40 antibody [MT86 + MT221] as coating antibody and GTX03013-02 IL12 p70 antibody [MT704] (Biotin) as detecting antibody.



For full product information, images and publications, please visit our [website](#).