

IL-12 p70 ELISA pair [C17.8/IL12-IV]

Cat. No. GTX03074

实验应用 ELISA, Sandwich ELISA

种属反应 Mouse

Package

1 pair

产品说明

Product Content

Cat No.	Type	宿主	偶联	Package
GTX02952-02	Detection antibody	Rat	Biotin	250 µg
GTX02954	Capture antibody	Rat	Unconjugated	250 µg

实验应用

产品说明

GTX02954 (C17.8) is able to detect native and recombinant mouse p40 subunit of IL-12 and IL-23. In combination with GTX02952-02 (IL12-IV), it recognizes IL-12 p70 heterodimer (p40 + p35).

属性

形式 Liquid


存放说明

Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.

注意事项

仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

产品图片

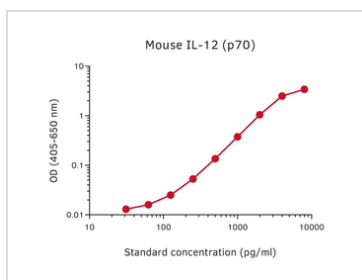

GTX03074 ELISA Image

Sandwich ELISA analysis of mouse IL-12 (p70) protein using GTX02954 IL12B / IL12 p40 antibody [C17.8] as coating antibody and GTX02952-02 IL12A / IL12 p35 antibody [IL12-IV] (Biotin) as detecting antibody.



For full product information, images and publications, please visit our [website](#).