

MUTYH antibody

Cat. No. GTX03359

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, ICC/IF, IHC-P
种属反应	Human, Mouse

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500-1:2000
ICC/IF	1:100-1:500
IHC-P	1:50-1:200

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 150mM NaCl, 50% Glycerol
保存剂	0.02% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

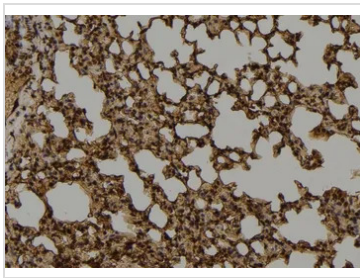


GTX03359 IHC-P Image

IHC-P analysis of mouse lung tissue section using GTX03359 MUTYH antibody.
Antigen retrieval : Heat mediated antigen retrieval step in citrate buffer was performed.
Dilution : 1:100

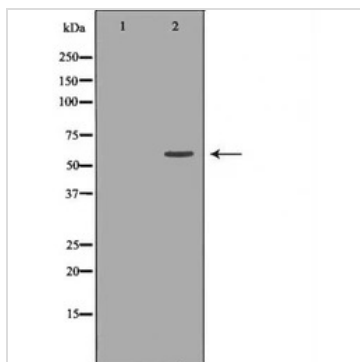

GTX03359 WB Image

WB analysis of HeLa whole cell lysate using GTX03359 MUTYH antibody. The lane on the left was treated with the antigen-specific peptide.

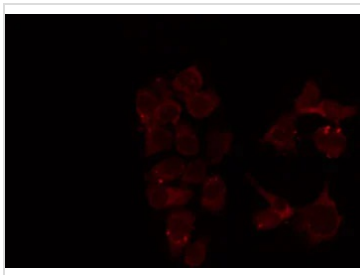

GTX03359 ICC/IF Image

ICC/IF analysis of PFA-fixed HeLa cells using GTX03359 MUTYH antibody.
Permeabilization : 0.1% Triton X-100
Dilution : 1:200



For full product information, images and publications, please visit our [website](#).