

AMFR antibody

Cat. No. GTX04315

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, ICC/IF, IHC-P, FCM
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	0.1-0.5 μ g/ml
ICC/IF	5 μ g/ml
IHC-P	0.5-1 μ g/ml
FCM	1-3 μ g/1x10 ⁶ cells

Note : Antigen retrieval by heat

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片

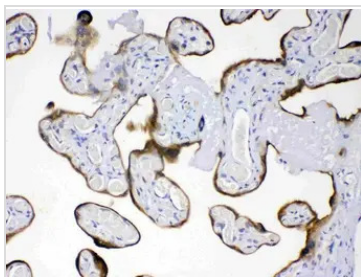
属性

形式	Liquid
存储溶液	0.2 mg Na ₂ HPO ₄ , 0.9 mg NaCl, 5 mg BSA
保存剂	0.05mg sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.5 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

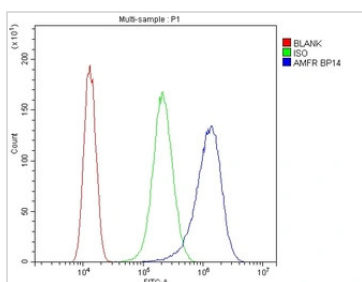
產品圖片

**GTx04315 IHC-P Image**

IHC-P analysis of human placenta tissue using GTx04315 AMFR antibody.

Antigen retrieval : Heat mediated antigen retrieval was performed in EDTA buffer (pH 8.0)

Dilution : 1 µg/ml

**GTx04315 FCM Image**

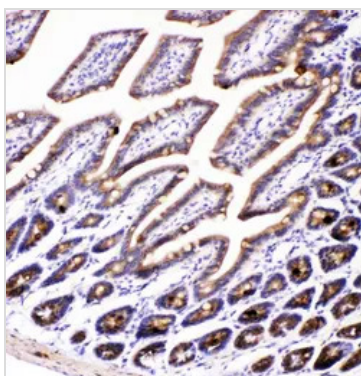
FACS analysis of SiHa cells using GTx04315 AMFR antibody.

Dilution : 1µg/1x10⁶ cells

Blue : Primary antibody

Green : Isotype control

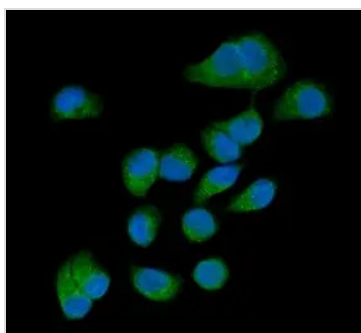
Red : Unlabelled sample

**GTx04315 IHC-P Image**

IHC-P analysis of mouse intestine tissue using GTx04315 AMFR antibody.

Antigen retrieval : Heat mediated antigen retrieval was performed in citrate buffer (pH 6.0) for 20 mins

Dilution : 1 µg/ml

**GTx04315 ICC/IF Image**

ICC/IF analysis of T47D cells using GTx04315 AMFR antibody.

Green : Primary antibody

Blue : DAPI

Antigen retrieval : IHC enzyme antigen retrieval reagent for 15 mins

Dilution : 5 µg/mL



For full product information, images and publications, please visit our [website](http://www.genetex.com).