

MDMX antibody

Cat. No. GTX04638

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P
种属反应	Human, Mouse

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	0.1-0.5 µg/mL
IHC-P	0.5-1 µg/mL

Note : By heat

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	0.2 mg Na ₂ HPO ₄ , 0.9 mg NaCl, 5 mg BSA
保存剂	0.05 mg Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.5 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

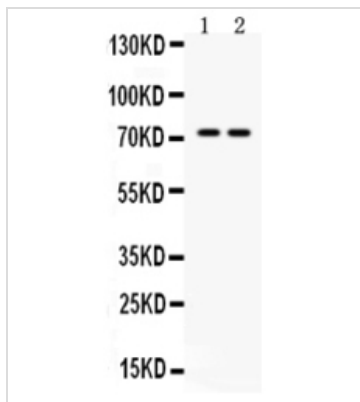

GTx04638 WB Image


WB analysis of various samples using GTx04638 MDMX antibody.

Lane 1 : Mouse testis tissue lysate

Lane 2 : 22RV1 whole cell lysate

Loading : 30 µg

Dilution : 0.5 µg/mL

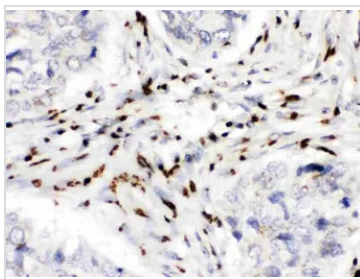

GTx04638 IHC-P Image

IHC-P analysis of human lung cancer tissue using GTx04638 MDMX antibody.

Antigen retrieval : Heat mediated EDTA buffer (pH 8.0)

Dilution : 1 µg/mL



For full product information, images and publications, please visit our [website](#).