

Pan phosphoserine/threonine antibody [ARC64485]

Cat. No. GTX05060

宿主	Rabbit
克隆	Monoclonal
同种型	IgG
实验应用	WB, IHC-P
种属反应	Species independent, Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500-1:1000
IHC-P	1:50-1:200

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	PBS, 0.05% BSA, 50% Glycerol
保存剂	0.09% Sodium Azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. Aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



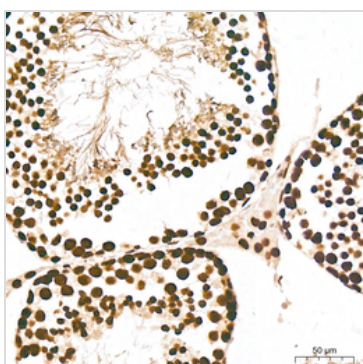
For full product information, images and publications, please visit our [website](#).

產品圖片

GTX05060 IHC-P Image

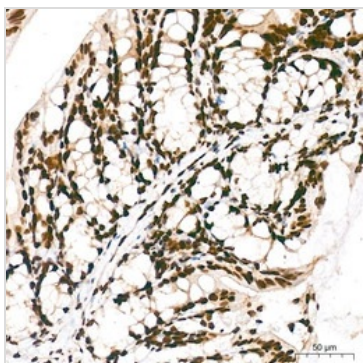
IHC-P analysis of human thyroid tissue using GTX05060 Pan phosphoserine/threonine antibody [ARC64485].

Antigen retrieval : High pressure mediated antigen retrieval performed with 0.01 M citrate buffer (pH 6.0)
Dilution : 1:200


GTX05060 IHC-P Image

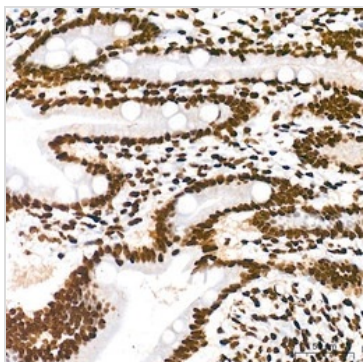
IHC-P analysis of mouse testis tissue using GTX05060 Pan phosphoserine/threonine antibody [ARC64485].

Antigen retrieval : High pressure mediated antigen retrieval performed with 0.01 M citrate buffer (pH 6.0)
Dilution : 1:200


GTX05060 IHC-P Image

IHC-P analysis of mouse colon tissue using GTX05060 Pan phosphoserine/threonine antibody [ARC64485].

Antigen retrieval : High pressure mediated antigen retrieval performed with 0.01 M citrate buffer (pH 6.0)
Dilution : 1:200


GTX05060 IHC-P Image

IHC-P analysis of human colon tissue using GTX05060 Pan phosphoserine/threonine antibody [ARC64485].

Antigen retrieval : High pressure mediated antigen retrieval performed with 0.01 M citrate buffer (pH 6.0)
Dilution : 1:200



For full product information, images and publications, please visit our [website](https://www.genetex.com).