

PAPP A antibody [10H9]

Cat. No. GTX10195

宿主	Mouse
克隆	Monoclonal
同种型	IgG2a
实验应用	IP, ELISA
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IP	Assay dependent
ELISA	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片

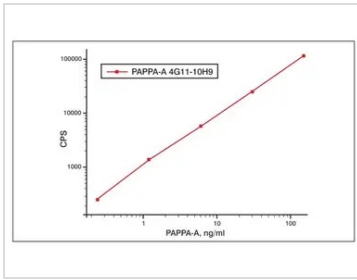
属性

形式	Liquid
存储溶液	PBS
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. Store at 4°C.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床，治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

**GTX10195 ELISA Image**

Calibration curve for PAPP-A sandwich immunoassay. Capture: PAPP A antibody [4G11] (GTX10192) biotinylated; Detection: PAPP A antibody [10H9] (GTX10195) labeled with stable Eu3+ -chelate; Antigen: htPAPP-A; Mixture of antibodies and antigen was incubated for 30 minutes at room temperature in streptavidin-coated plates.



For full product information, images and publications, please visit our [website](#).