

Destrin antibody

Cat. No. GTX11072

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, ICC/IF, IHC-P
种属反应	Human, Mouse, Rat, Dog

引用文献 (2)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:1500
ICC/IF	1:100
IHC-P	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	PBS, 1% BSA
保存剂	15mM Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	~0.5-1.0 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_381235
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。


For full product information, images and publications, please visit our [website](#).

產品圖片

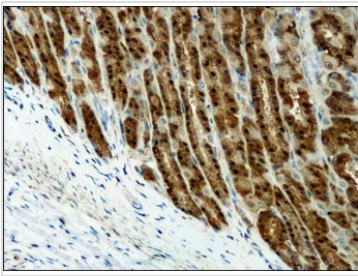


GTX11072 IHC-P Image

IHC-P analysis of rat stomach tissue using GTX11072 Destrin antibody at 0.02 µg/mL.

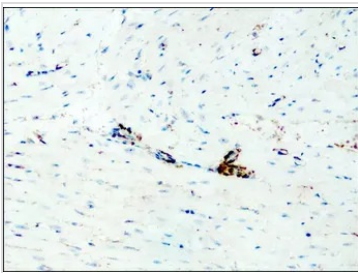

GTX11072 IHC-P Image

IHC-P analysis of mouse heart tissue using GTX11072 Destrin antibody at 0.02 µg/mL.

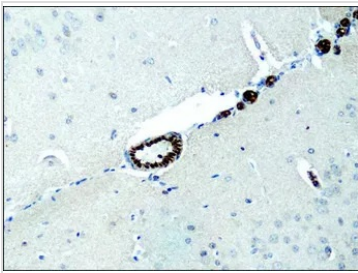

GTX11072 IHC-P Image

IHC-P analysis of mouse brain tissue using GTX11072 Destrin antibody at 0.02 µg/mL.



For full product information, images and publications, please visit our [website](#).