

gamma Tubulin antibody [GTU-88]

Cat. No. GTX11316

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	WB, ICC/IF, ELISA
种属反应	Human, Mouse, Rat, Drosophila, Dog, Hamster, Chicken, Monkey, Xenopus

引用文献 (5)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:10,000
ICC/IF	1:5,000-1:10,000
ELISA	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	Ascites
保存剂	15mM Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
偶联	Unconjugated
RRID	AB_384215
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。


For full product information, images and publications, please visit our [website](#).

產品圖片

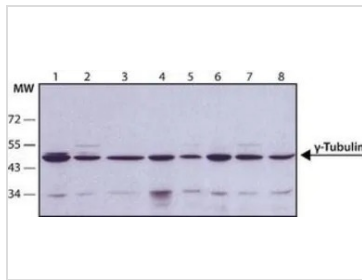


GTX11316 WB Image

WB analysis of (1) HeLa (2) U87 (3) COS7 (4) P19 (5) rat2 (6) CHO (7) MDBK (8) MDCK cell lysate using GTX11316 gamma Tubulin antibody [GTU-88] at 1:10,000.

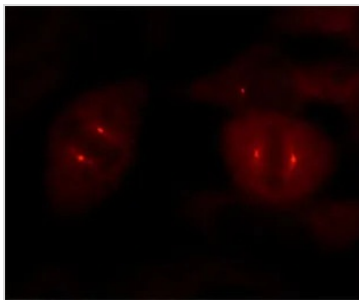

GTX11316 ICC/IF Image

ICC/IF analysis of HeLa cells using GTX11316 gamma Tubulin antibody [GTU-88] at 1:5,000. Cells were fixed and permeabilized with methanol followed by acetone.

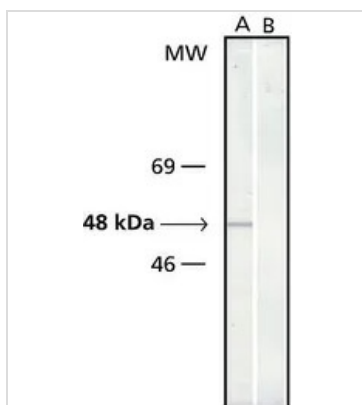


GTX11316 WB Image

WB analysis of chicken fibroblasts using GTX11316 gamma Tubulin antibody [GTU-88].

Lane A : Antibody dilution 1:20,000

Lane B : Negative control (secondary antibody only)

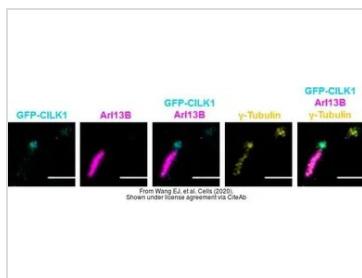

GTX11316 ICC/IF Image

The data was published in the journal Cells in 2020. [PMID: 32178256](https://pubmed.ncbi.nlm.nih.gov/32178256/)



For full product information, images and publications, please visit our [website](https://www.genetex.com).