

Avian Influenza A virus H5N3 HA (Hemagglutinin) antibody

Cat. No. GTX127299

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB
种属反应	Avian Influenza A virus

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500-1:3000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 20% Glycerol
保存剂	0.01% Thimerosal
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_2885634
注意事项	仅供实验室使用。不适用于人类或动物的任何临床，治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

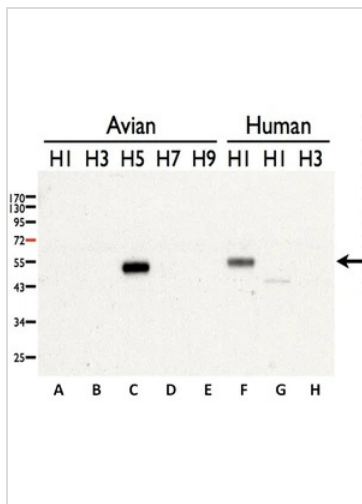

GTX127299 WB Image

Avian Influenza A Virus H5N3 Hemagglutinin (HA) antibody detects H5 (Influenza A) protein by Western blot analysis.

- A. MDCK cell lysates infected with 4 HAU/ml of A/duck/Alberta/35/1976-like, H1N1
- B. MDCK cell lysates infected with 4 HAU/ml of A/duck/Ukraine/1963-like, H3N8
- C. MDCK cell lysates infected with 4 HAU/ml of A/mallard/Miyagi/53/1976-like, H5N3
- D. MDCK cell lysates infected with 4 HAU/ml of A/duck/Taiwan/33/1993-like, H7N7
- E. MDCK cell lysates infected with 4 HAU/ml of A/turkey/Wisconsin/1/1966-like, H9N2
- F. MDCK cell lysates infected with 4 HAU/ml of A/PR/8/34, H1N1
- G. MDCK cell lysates infected with 4 HAU/ml of A/Taiwan/126/2009, pdmH1N1
- H. MDCK cell lysates infected with 4 HAU/ml of A/Taiwan/3446/2002, H3N2

12 % SDS-PAGE

Avian Influenza A Virus H5N3 Hemagglutinin (HA) antibody (GTX127299) dilution: 1:1000



For full product information, images and publications, please visit our [website](http://www.genetex.com).