

Connexin 40 antibody

Cat. No. GTX128055

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB
种属反应	Human, Mouse

引用文献 (2)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500-1:3000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片

属性


形式	Liquid
存储溶液	PBS, 20% Glycerol
保存剂	0.01% Thimerosal
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_2885692
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

Date 2025 / 05 / 10 Page 1 of 2

產品圖片

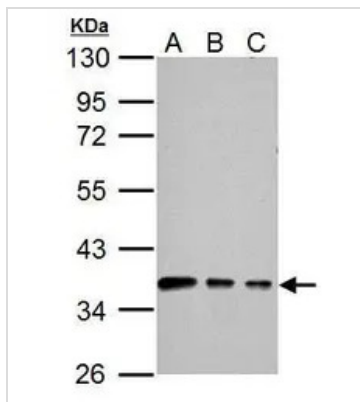

GTX128055 WB Image

Connexin 40 antibody detects Connexin 40 protein by Western blot analysis.

A. 30 μ g 293T whole cell lysate/extract

B. 30 μ g A431 whole cell lysate/extract

C. 30 μ g HepG2 whole cell lysate/extract

10 % SDS-PAGE

Connexin 40 antibody (GTX128055) dilution: 1:1000



For full product information, images and publications, please visit our [website](#).