

GPR81 antibody

Cat. No. GTX12960

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	IHC-P
种属反应	Human, Rat, Rabbit, Dog, Hamster, Monkey, Bat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IHC-P	3 µg/ml

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS
保存剂	0.1% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

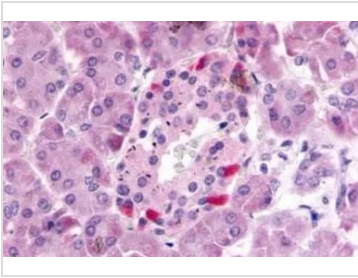


GTX12960 IHC-P Image

IHC-P analysis of human pancreas islet of langerhans tissue using GTX12960 GPR81 antibody.
Antigen retrieval : Heat-induced antigen retrieval

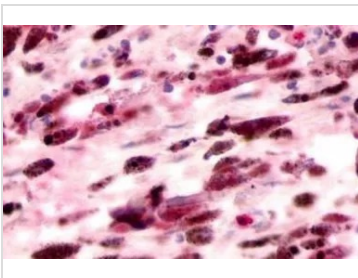


GTX12960 IHC-P Image

IHC-P analysis of human skin, melanoma tissue using GTX12960 GPR81 antibody.
Antigen retrieval : Heat-induced antigen retrieval

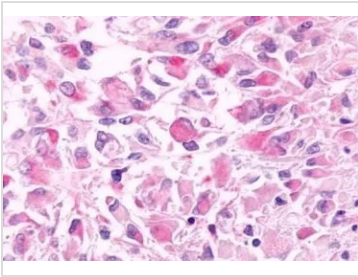


GTX12960 IHC-P Image

IHC-P analysis of brain, glioblastoma tissue using GTX12960 GPR81 antibody.

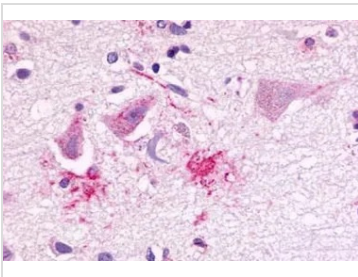

GTX12960 IHC-P Image

IHC-P analysis of brain, alzheimer's, neurons and senile plaque tissue using GTX12960 GPR81 antibody.

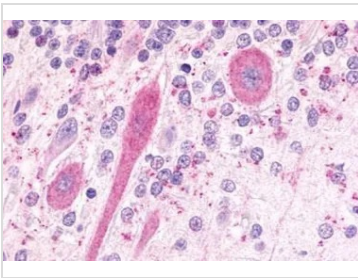

GTX12960 IHC-P Image

IHC-P analysis of brain, cerebellum, purkinje cell layer tissue using GTX12960 GPR81 antibody.



For full product information, images and publications, please visit our [website](https://www.genetex.com).