

MAS1 antibody

Cat. No. GTX13032

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	IHC-P
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IHC-P	15 µg/ml

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片

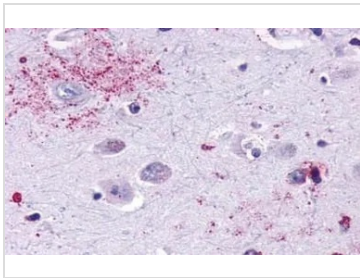
属性

形式	Liquid
存储溶液	PBS
保存剂	0.1% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。

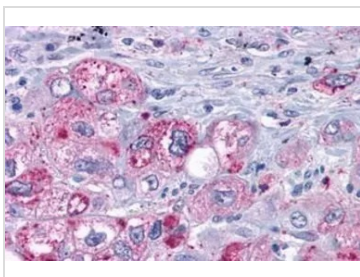


For full product information, images and publications, please visit our [website](#).

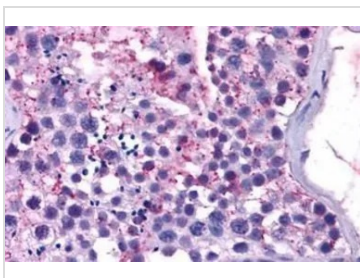
產品图片


GTX13032 IHC-P Image

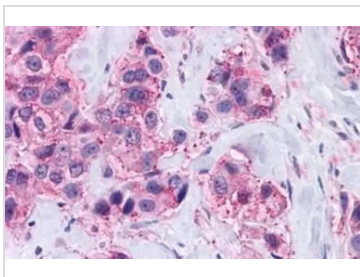
IHC-P analysis of human brain, neurons and glia tissue using GTX13032 MAS1 antibody.
Antigen retrieval : Heat-induced antigen retrieval


GTX13032 IHC-P Image

IHC-P analysis of human lung, adenocarcinoma tissue using GTX13032 MAS1 antibody.
Antigen retrieval : Heat-induced antigen retrieval


GTX13032 IHC-P Image

IHC-P analysis of human testis tissue using GTX13032 MAS1 antibody.
Antigen retrieval : Heat-induced antigen retrieval


GTX13032 IHC-P Image

IHC-P analysis of human breast, carcinoma tissue using GTX13032 MAS1 antibody.
Antigen retrieval : Heat-induced antigen retrieval



For full product information, images and publications, please visit our [website](#).