

Hsp90 antibody [AC88]

Cat. No. GTX13492

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	WB, ICC/IF, IHC-P, FCM, IP
种属反应	Human, Mouse, Rat, Rabbit, Sheep, Bovine, Dog, Hamster, Guinea pig, Pig, Monkey, Carp, Rainbow trout, Fish

引用文献 (1)

实验应用


应用说明

IHC-P: Use at an assay dependent dilution. When used in immunoprecipitation, this antibody does not appear to react well with Hsp90 that is bound to steroid receptors or to pp60v-src of Rouse sarcoma virus. WB: Use at a concentration of 1 µg/ml. Detects a band of approximately 90 kDa (predicted molecular weight: 84.7 kDa(alpha) and 83.2 kDa(beta)). Optimal dilutions/concentrations should be determined by the end user.

属性

形式	Liquid
存储溶液	PBS, 50% Glycerol, 0.1mM PMSF
保存剂	No preservative
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

产品图片

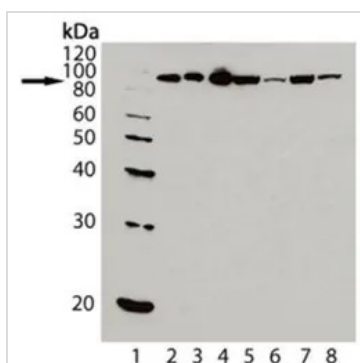


GTX13492 WB Image

Western blot analysis of HSP90, mAb (AC88) : Lane 1: MW Marker; Lane 2: HSP90 (human), (native) ; Lane 3: HSP90beta, (human) (recombinant) ; Lane 4: HSP90alpha (human), (recombinant) ; Lane 5: HeLa (heat shocked); Lane 6: 3T3 (heat shocked) ; Lane 7: PC-12 (heat shocked) ; Lane 8: CHO-K1 (heat shocked).



For full product information, images and publications, please visit our [website](#).

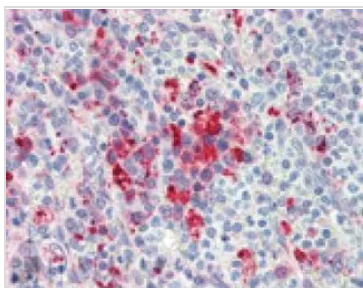
GTX13492 IHC-P Image

Immunohistochemistry(paraffin sections) analysis of human spleen tissue stained with HSP90, mAb (AC88) at 10µg/ml.



For full product information, images and publications, please visit our [website](#).