

A2B5 antibody [105]

Cat. No. GTX14538

宿主	Mouse
克隆	Monoclonal
同种型	IgM
实验应用	ICC/IF
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
ICC/IF	5-10 µg/mL

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 5% Trehalose
保存剂	No preservative
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

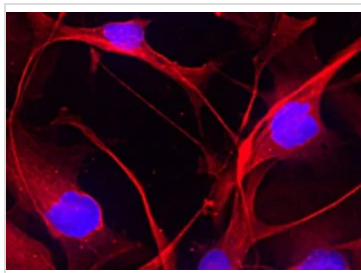

GTX14538 ICC/IF Image

ICC/IF analysis of A172 Cells using GTX14538 A2B5 antibody [105].

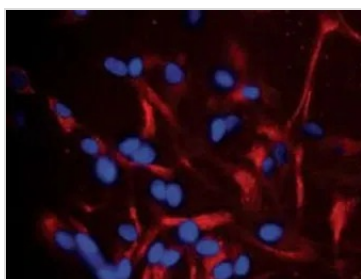

GTX14538 ICC/IF Image

ICC/IF analysis of human neural progenitors cells using GTX14538 A2B5 antibody [105].

Red : Primary antibody

Blue : DAPI



For full product information, images and publications, please visit our [website](#).