

CD44 antibody [IM7] (Biotin)

Cat. No. GTX19629

宿主	Rat
克隆	Monoclonal
同种型	IgG2b
实验应用	IHC-Fr, FCM, IP
种属反应	Human, Bovine, Cat, Dog, Pig, Horse

引用文献 (4)

实验应用


应用说明

FACS: Use 1 μ g for 10⁶ cells. IHC-Fr (acetone fixed): Use at an assay dependent dilution. IP: Use at an assay dependent dilution. Optimal dilutions/concentrations should be determined by the end user.

属性

形式	Liquid
存储溶液	PBS, 1% BSA
保存剂	0.1% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. Store at 4°C. DO NOT FREEZE.
浓度	0.5 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Biotin
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。

产品图片

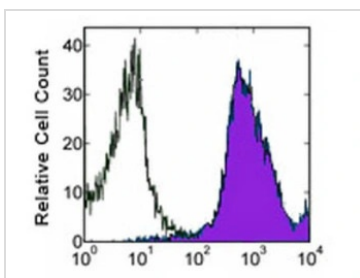

GTX19629 FCM Image

Staining of C57 splenocytes with Anti-Mouse CD44: Biotin (purple histogram) or Rat IgG2b Isotype Control: Biotin (open histogram) followed by UltraAvidin: PE



For full product information, images and publications, please visit our [website](#).

Date 2026 / 06 / 06 Page 1 of 1