

CCR8 antibody

Cat. No. GTX21666

宿主	Goat
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P, FCM, ELISA, Neutralizing/Inhibition
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500
IHC-P	1:125
FCM	2-6 μ l/10 ⁶ cells/50-100 μ l
ELISA	1:50,000
Neutralizing/Inhibition	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	10mM KHPO ₄ , 140mM NaCl, 1mg/ml BSA
保存剂	0.1% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

For full product information, images and publications, please visit our [website](#).

產品圖片

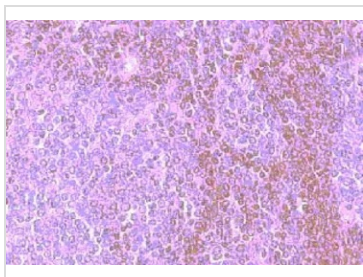

GTX21666 IHC-P Image

IHC-P analysis of human spleen tissue using GTX21666 CCR8 antibody.



For full product information, images and publications, please visit our [website](#).