

Androgen Receptor antibody [AN1-15]

Cat. No. GTX22742

宿主	Rat
克隆	Monoclonal
同种型	IgG
实验应用	WB, ICC/IF, IHC-P, IHC-Fr, IP, Gel supershift assays, IHC
种属反应	Human, Mouse, Rat, Primate

引用文献 (4)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:50-1:500
ICC/IF	5 µg/ml
IHC-P	2-4 µg/mL
IHC-Fr	2-4 µg/mL
IP	Assay dependent
Gel supershift assays	Assay dependent
IHC	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	PBS
保存剂	0.05% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_367517
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

For full product information, images and publications, please visit our [website](#).

產品圖片

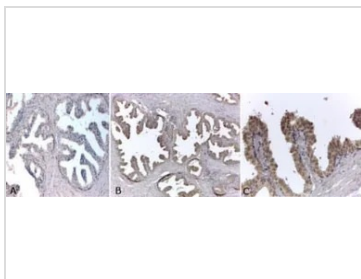

GTX22742 IHC-P Image

IHC-P analysis of human prostate tissue using GTX22742 Androgen Receptor antibody [AN1-15].


A : isotype control

B : Primary antibody

C : Primary antibody with magnified section


Antigen retrieval : HEIR with a buffer (pH 6.2)

Dilution : 2 µg/ml

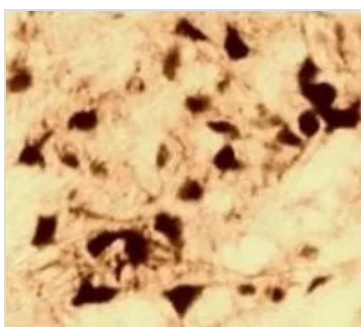

GTX22742 IHC-P Image

IHC-P analysis of primate brain tissue using GTX22742 Androgen Receptor antibody [AN1-15].

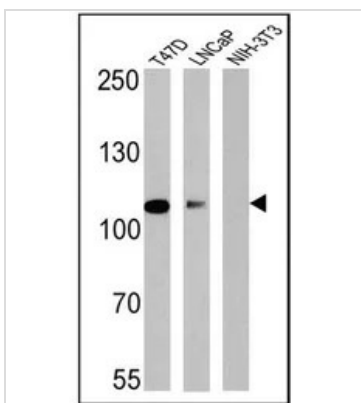

GTX22742 WB Image

WB analysis of 25 µg of T47D (lane 1), LNCaP (lane 2) and NIH-3T3 (lane 3) lysates using GTX22742 Androgen Receptor antibody [AN1-15].

Dilution : 1:100



For full product information, images and publications, please visit our [website](https://www.genetex.com).