

Glucocorticoid Receptor antibody [BuGR2]

Cat. No. GTX22768

宿主	Mouse
克隆	Monoclonal
同种型	IgG2a
实验应用	WB, ICC/IF, IHC-P, FCM, IP, ELISA, ChIP assay, Gel supershift assays
种属反应	Human, Mouse, Rat, Yeast, Rabbit, Sheep, Guinea pig, Xenopus, Primate

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	5 µg/ml
ICC/IF	1:50-1:500
IHC-P	5 µg/ml
FCM	1-2 µg/test
IP	Assay dependent
ELISA	Assay dependent
ChIP assay	Assay dependent
Gel supershift assays	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS
保存剂	0.05% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

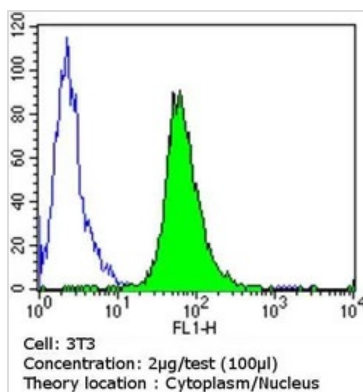


GTX22768 FCM Image

FACS analysis of NIH/3T3 cells using GTX22768 Glucocorticoid Receptor antibody [BuGR2] compared to an isotype control (blue).

Dilution : 2 ug/test for 60 min at room temperature


Fixation : 2% paraformaldehyde

Permeabilization : 90% methanol and incubated for 10 minutes at room temperature

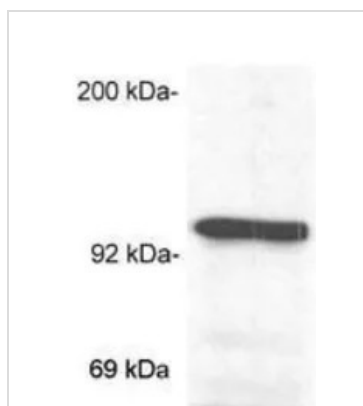

GTX22768 WB Image

WB analysis of mouse liver extract using GTX22768 Glucocorticoid Receptor antibody [BuGR2].

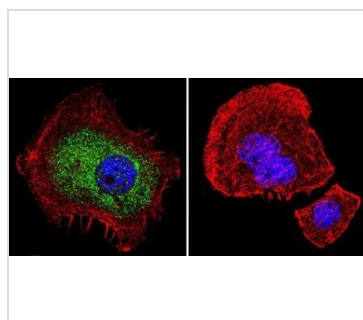

GTX22768 ICC/IF Image

ICC/IF analysis of A549 cells using GTX22768 Glucocorticoid Receptor antibody [BuGR2]. Cells were probed without (right) or with(left) an antibody.

Green : Primary antibody

Blue : Nuclei

Red : Actin

Fixation : formaldehyde

Dilution : 1:100 overnight at 4 °C



For full product information, images and publications, please visit our [website](https://www.genetex.com).