

Goat Anti-Rabbit IgG (F(ab')₂) antibody, F(ab')₂ fragment, pre-adsorbed (FITC)

Cat. No. GTX26108

宿主	Goat
克隆	Polyclonal
同种型	IgG F(ab') ₂
实验应用	ICC/IF, FCM, Dot, ELISA
种属反应	Rabbit

Package
250 µg

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
ICC/IF	1:1000-1:5000
FCM	1:500-1:2500
Dot	Assay dependent
ELISA	1:10000-1:50000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片

产品说明

Pre-adsorbed with Bovine, Horse, Human, Mouse, Rat and Sheep serum proteins. May react with immunoglobulins from other species.

属性

形式	Liquid
存储溶液	20mM Potassium Phosphate, 150mM NaCl, 1% BSA
保存剂	0.01% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. Protect from light.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
免疫原	Rabbit IgG F(ab') ₂ fragment
纯化方式	Purified by antigen-affinity chromatography using Rabbit IgG coupled to agarose beads followed by solid phase adsorptions to remove any unwanted reactivities, pepsin digestion and chromatographic separation. From serum


For full product information, images and publications, please visit our [website](#).

Date 2025 / 05 / 03 Page 1 of 2

偶联 Fluorescein isothiocyanate (FITC) [Wavelength](#)

注意事项 仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。

產品图片

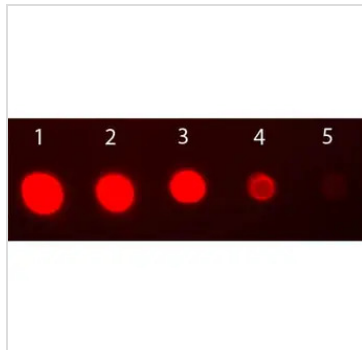

GTx26108 Dot Image

Dot blot analysis of rabbit IgG using GTx26108 Goat Anti-Rabbit IgG (F(ab')₂) antibody, F(ab')₂ fragment, pre-adsorbed (FITC).

Lane 1 : 50 ng

Lane 2 : 16.67 ng

Lane 3 : 5.56 ng

Lane 4 : 1.85 ng

Lane 5 : 0.62 ng

Dilution : 1:1000



For full product information, images and publications, please visit our [website](#).