

EGF antibody

Cat. No. GTX29695

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P, ELISA, Neutralizing/Inhibition, Sandwich ELISA
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	0.1 - 0.2 µg/mL
IHC-P	Assay dependent
ELISA	0.1 - 0.2 µg/mL
Neutralizing/Inhibition	0.25 - 0.35 µg/mL
Sandwich ELISA	0.5 - 2.0 ug/ml

Note : To yield one-half maximal inhibition [ND₅₀] of the biological activity of hEGF (2.0 ng/mL)
Capture : GTX29695, Detection : Biotinylated Anti-Human EGF antibody

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS
保存剂	No preservative
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

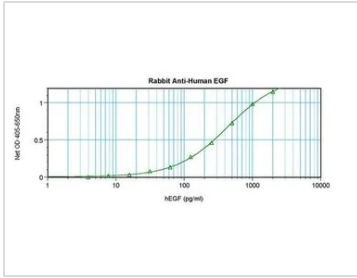


GTX29695 ELISA Image

ELISA analysis of human EGF recombinant protein(0.2 - 0.4 ng/well) using GTX29695 EGF antibody(capture antibody) at 0.5 - 2.0 ug/ml and Biotinylated Anti-Human EGF as a detection antibody.

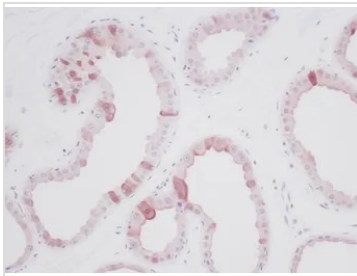


GTX29695 IHC-P Image

IHC-P analysis of human skin tissue using GTX29695 EGF antibody.
Working concentration : 1.0 ug/ml - 2.5 ug/ml
Antigen retrieval : proteinase K.

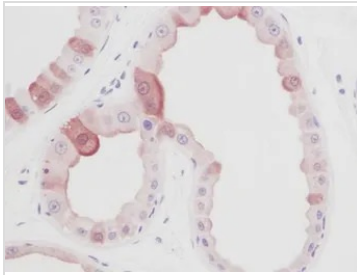

GTX29695 IHC-P Image

IHC-P analysis of human skin tissue using GTX29695 EGF antibody.
Working concentration : 1.0 ug/ml - 2.5 ug/ml
Antigen retrieval : proteinase K.

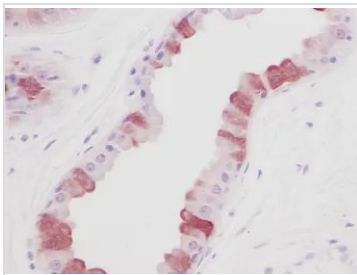

GTX29695 IHC-P Image

IHC-P analysis of human skin tissue using GTX29695 EGF antibody.
Working concentration : 1.0 ug/ml - 2.5 ug/ml
Antigen retrieval : proteinase K.



For full product information, images and publications, please visit our [website](https://www.genetex.com).