

GTX300067

TLRs Antibody Panel (intracellular menbrane)

Product Content

Cat No	Product Name	种属反应	实验应用	Package
GTX100726	TLR9 antibody [N2], N-term	Human, Mouse	WB, IHC-P	25 µl
GTX125910	TLR7 antibody	Human	WB, IHC-P	25 µl
GTX213110-01	Goat Anti-Rabbit IgG antibody (HRP)	Rabbit	WB, IHC-P, ELISA	25 µl
GTX638477	TLR3 antibody [HL2314]	Human, Mouse, Rat, Cat, Dog	WB, IHC-P	25 µl

注意事项

For In vitro laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.



For full product information, images and publications, please visit our [website](#).

TLR9 antibody [N2], N-term

Cat. No. GTX100726

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P
种属反应	Human, Mouse

引用文献 (1)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500-1:3000
IHC-P	1:100-1:1000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	0.1M Tris, 0.1M Glycine, 10% Glycerol
保存剂	0.01% Thimerosal
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.4 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_2038134
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

For full product information, images and publications, please visit our [website](#).

產品圖片

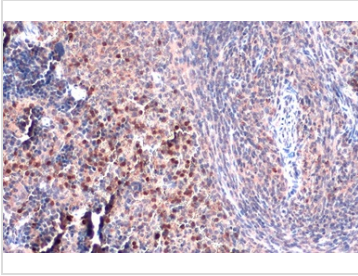


GTX100726 IHC-P Image

TLR9 antibody [N2], N-term detects TLR9 protein by immunohistochemical analysis.

Sample: Paraffin-embedded mouse spleen.

TLR9 stained by TLR9 antibody [N2], N-term (GTX100726) diluted at 1:500.


Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

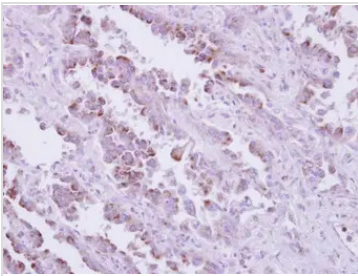

GTX100726 IHC-P Image

Immunohistochemical analysis of paraffin-embedded human lung cancer, using TLR9(GTX100726) antibody at 1:250 dilution.

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

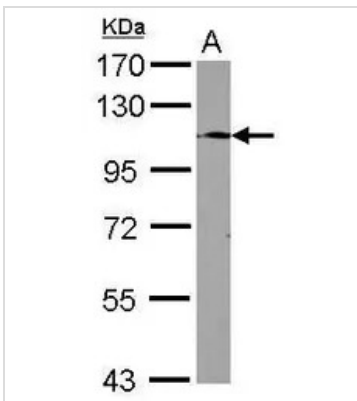


GTX100726 WB Image

Sample (30 ug of whole cell lysate)

A: THP-1

7.5% SDS PAGE

GTX100726 diluted at 1:1000

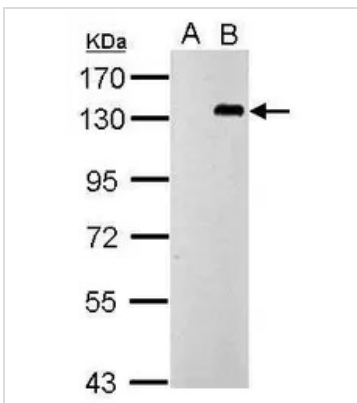

GTX100726 WB Image

Sample (30 ug of whole cell lysate)

A: Non-transfected 293T lysates

B: TLR9 transfected 293T lysates

7.5% SDS PAGE

GTX100726 diluted at 1:2000



For full product information, images and publications, please visit our [website](http://www.genetex.com).

TLR7 antibody

Cat. No. GTX125910

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P
种属反应	Human

引用文献 (1)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:1000-1:10000
IHC-P	1:100-1:1000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 20% Glycerol
保存剂	0.01% Thimerosal
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_11176526
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

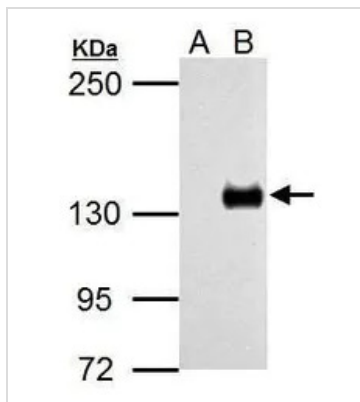

GTx125910 WB Image


Sample (30 ug of whole cell lysate)

A: Non-transfected 293T lysates

B: TLR7 transfected 293T lysates

5% SDS PAGE

GTx125910 diluted at 1:5000

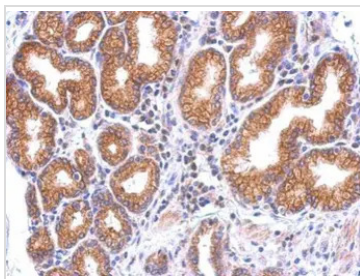

GTx125910 IHC-P Image

Immunohistochemical analysis of paraffin-embedded human colon carcinoma, using TLR7 antibody (GTx125910) antibody at 1:500 dilution.

Antigen Retrieval: Trilogy™ (EDTA based, pH 8.0) buffer, 15min



For full product information, images and publications, please visit our [website](https://www.genetex.com).

Goat Anti-Rabbit IgG antibody (HRP)

Cat. No. GTX213110-01

宿主	Goat
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P, ELISA
种属反应	Rabbit

引用文献 (562)
Package
1 ml

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	Assay dependent
IHC-P	1:100-1:1000
ELISA	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	0.05M Tris, 0.15M NaCl, 1%BSA
保存剂	0.025% ProClin 300
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.53 mg/ml (Please refer to the vial label for the specific concentration.)
免疫原	Highly purified whole rabbit IgG
纯化方式	Purified by antigen-affinity chromatography.
偶联	Horseradish peroxidase(HRP)
RRID	AB_10618573
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

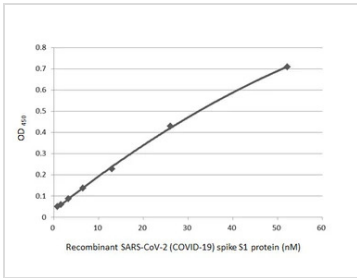**GTX213110-01 ELISA Image**

Indirect ELISA analysis performed by coating plate with recombinant SARS-CoV-2 (COVID-19) Spike S1 protein, His tag (active) (GTX135817-pro) (52.15-0.81 nM). Coated protein was probed with SARS-CoV-2 (COVID-19) Spike S1 antibody [HL134] (GTX635671) (1 µg/mL). Goat anti-rabbit IgG antibody (HRP) (GTX213110-01) (1:10000) was used to detect bound primary antibody.

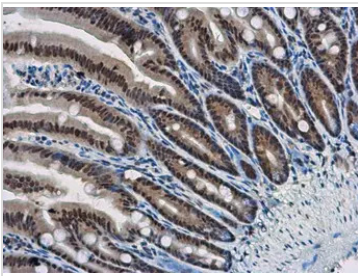
**GTX213110-01 IHC-P Image**

WBP11 antibody detects WBP11 protein at nucleus in mouse intestine by immunohistochemical analysis. Sample: Paraffin-embedded mouse intestine.


WBP11 antibody (GTX118654) diluted at 1:500.

The signal was developed by Rabbit IgG antibody (HRP) (GTX213110-01)

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

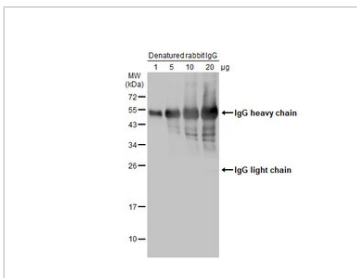**GTX213110-01 WB Image**

Various amounts of denatured rabbit IgG protein were separated by 12% SDS-PAGE, and the membrane was blotted with HRP-conjugated anti-rabbit IgG antibody (GTX213110-01) diluted at 1:5000.

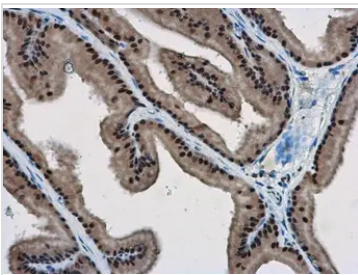**GTX213110-01 IHC-P Image**


WBP11 antibody detects WBP11 protein at nucleus in rat prostate by immunohistochemical analysis.

Sample: Paraffin-embedded rat prostate.

WBP11 antibody (GTX118654) diluted at 1:500.

The signal was developed by Rabbit IgG antibody (HRP) (GTX213110-01).


Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

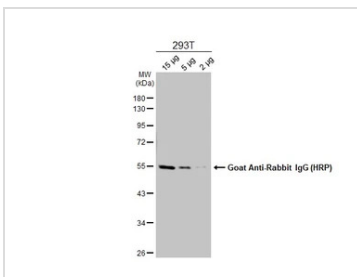**GTX213110-01 WB Image**

Various whole cell extracts were separated by 10% SDS-PAGE, and the membrane was blotted with alpha Tubulin antibody (GTX112141) diluted at 1:10000. The HRP-conjugated anti-rabbit IgG antibody (GTX213110-01) was used to detect the primary antibody.




For full product information, images and publications, please visit our [website](https://www.genetex.com).

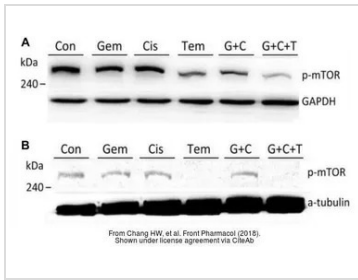


GTx213110-01 WB Image

The data was published in the journal Front Pharmacol in 2018. [PMID: 30087612](https://pubmed.ncbi.nlm.nih.gov/30087612/)

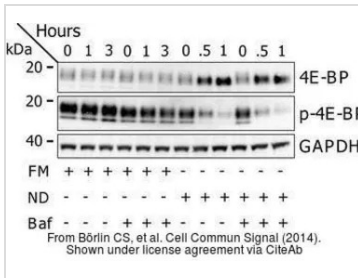


GTx213110-01 WB Image

The data was published in the journal Cell Commun Signal in 2014. [PMID: 25214434](https://pubmed.ncbi.nlm.nih.gov/25214434/)

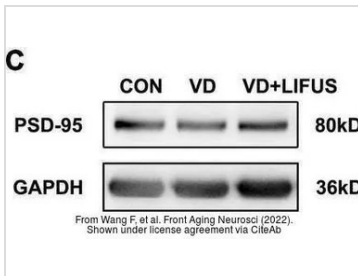


GTx213110-01 WB Image

The data was published in the 2022 in Front Aging Neurosci. [PMID: 35264943](https://pubmed.ncbi.nlm.nih.gov/35264943/)

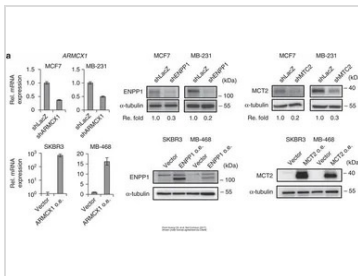

GTx213110-01 WB Image

The data was published in the 2017 in Nat Commun. [PMID: 28281525](https://pubmed.ncbi.nlm.nih.gov/28281525/)

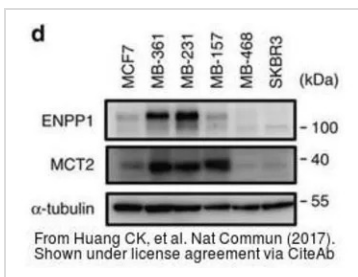


GTx213110-01 WB Image

The data was published in the 2017 in Nat Commun. [PMID: 28281525](https://pubmed.ncbi.nlm.nih.gov/28281525/)




For full product information, images and publications, please visit our [website](https://www.genetex.com).

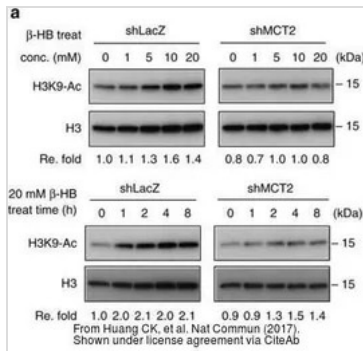


GTx213110-01 WB Image

The data was published in the 2017 in Nat Commun. [PMID: 28281525](#)

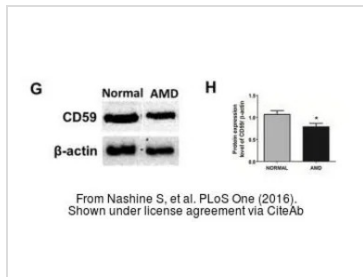


GTx213110-01 WB Image

The data was published in the 2016 in PLoS One. [PMID: 27486856](#)

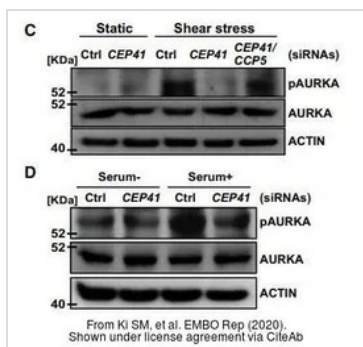

GTx213110-01 WB Image

The data was published in the journal EMBO Rep in 2020. [PMID: 31885126](#)

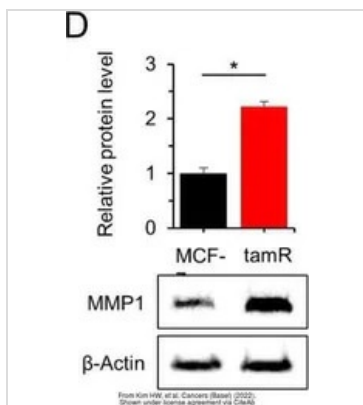


GTx213110-01 WB Image

The data was published in the 2022 in Cancers (Basel). [PMID: 35267540](#)




For full product information, images and publications, please visit our [website](#).

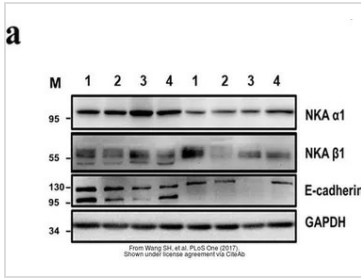


GTx213110-01 WB Image

The data was published in the journal PLoS One in 2017. [PMID: 28832634](https://pubmed.ncbi.nlm.nih.gov/28832634/)

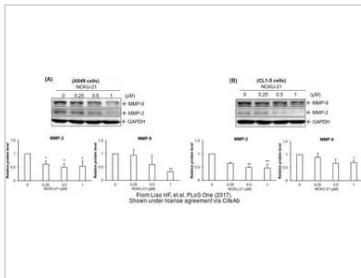


GTx213110-01 WB Image

The data was published in the journal PLoS One in 2017. [PMID: 28945763](https://pubmed.ncbi.nlm.nih.gov/28945763/)

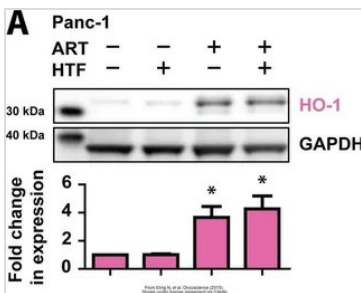


GTx213110-01 WB Image

The data was published in the journal Oncoscience in 2015. [PMID: 26097885](https://pubmed.ncbi.nlm.nih.gov/26097885/)

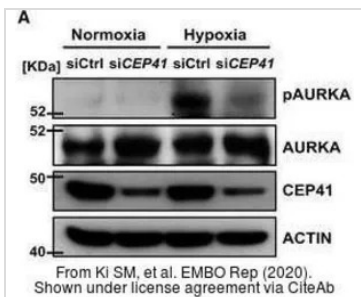

GTx213110-01 WB Image

The data was published in the journal EMBO Rep in 2020. [PMID: 31885126](https://pubmed.ncbi.nlm.nih.gov/31885126/)

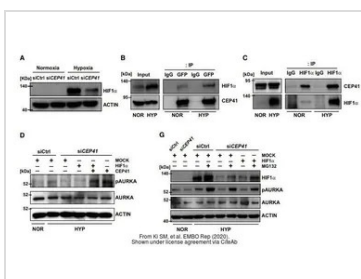

GTx213110-01 WB Image

The data was published in the journal EMBO Rep in 2020. [PMID: 31885126](https://pubmed.ncbi.nlm.nih.gov/31885126/)



For full product information, images and publications, please visit our [website](https://www.genetex.com).

TLR3 antibody [HL2314]

Cat. No. GTX638477

宿主	Rabbit
克隆	Monoclonal
同种型	IgG
实验应用	WB, IHC-P
种属反应	Human, Mouse, Rat, Cat, Dog

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500-1:3000
IHC-P	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS
保存剂	No preservative
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

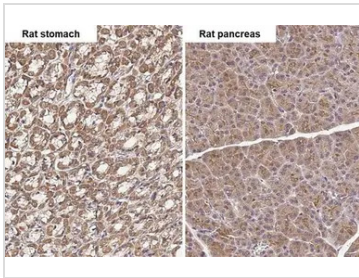


GTX638477 IHC-P Image

TLR3 antibody [HL2314] detects TLR3 protein by immunohistochemical analysis.

Sample: Paraffin-embedded rat tissues.

TLR3 stained by TLR3 antibody [HL2314] (GTX638477) diluted at 1:100.

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

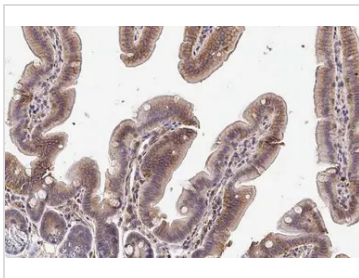


GTX638477 IHC-P Image

TLR3 antibody [HL2314] detects TLR3 protein at cytoplasm by immunohistochemical analysis.

Sample: Paraffin-embedded mouse intestine.

TLR3 stained by TLR3 antibody [HL2314] (GTX638477) diluted at 1:100.

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

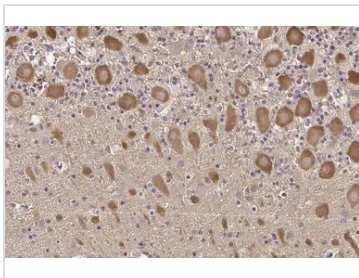


GTX638477 IHC-P Image

TLR3 antibody [HL2314] detects TLR3 protein at cytoplasm by immunohistochemical analysis.

Sample: Paraffin-embedded dog brain.


TLR3 stained by TLR3 antibody [HL2314] (GTX638477) diluted at 1:100.

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

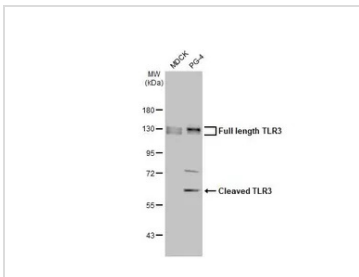

GTX638477 WB Image

Various whole cell extracts (30 µg) were separated by 7.5% SDS-PAGE, and the membrane was blotted with TLR3 antibody [HL2314] (GTX638477) diluted at 1:1000. The HRP-conjugated anti-rabbit IgG antibody (GTX213110-01) was used to detect the primary antibody.

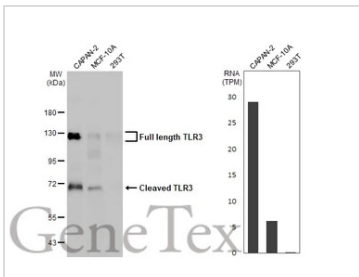

GTX638477 WB Image

Various whole cell extracts (30 µg) were separated by 7.5% SDS-PAGE, and the membrane was blotted with TLR3 antibody [HL2314] (GTX638477) diluted at 1:1000. The HRP-conjugated anti-rabbit IgG antibody (GTX213110-01) was used to detect the primary antibody, and the signal was developed with Trident femto Western HRP Substrate. Corresponding RNA expression data for the same cell lines are based on Human Protein Atlas program.



For full product information, images and publications, please visit our [website](https://www.genetex.com).