

C16orf80 antibody

Cat. No. GTX32483

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, ICC/IF
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:500 - 1:2000
ICC/IF	1:50 - 1:100

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 50% Glycerol
保存剂	0.05% ProClin 300
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床，治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品图片

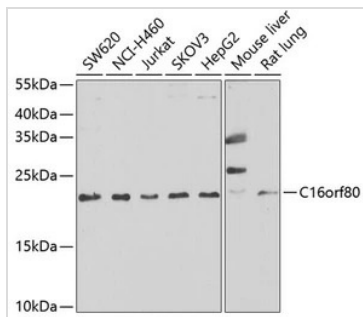

GTX32483 WB Image

WB analysis of various sample lysates using GTX32483 C16orf80 antibody. The signal was developed with ECL plus-Enhanced.
 Dilution : 1:1000
 Loading : 25µg per lane

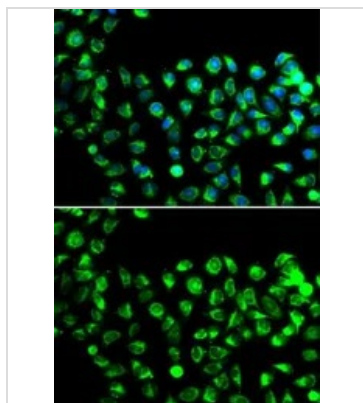

GTX32483 ICC/IF Image

ICC/IF analysis of A549 cells using GTX32483 C16orf80 antibody.
 Blue : DAPI



For full product information, images and publications, please visit our [website](#).