

## MTF2 antibody

Cat. No. GTX33330

|      |                   |
|------|-------------------|
| 宿主   | Rabbit            |
| 克隆   | Polyclonal        |
| 同种型  | IgG               |
| 实验应用 | WB, ICC/IF        |
| 种属反应 | Human, Mouse, Rat |

## 实验应用

## 应用说明

\*最佳稀释倍数与浓度应由研究人员确认

| Suggested dilution | Recommended dilution |
|--------------------|----------------------|
| WB                 | 1:500 - 1:2000       |
| ICC/IF             | 1:50 - 1:100         |

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

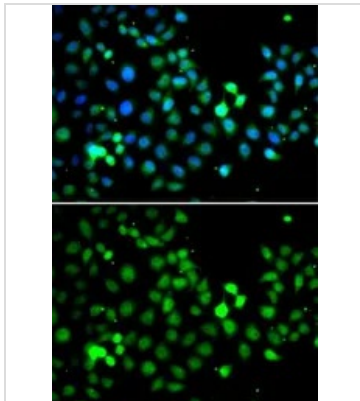
IHC-Fr: 冰冻切片

## 属性

|      |                                                                                                                                                                                                                            |
|------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 形式   | Liquid                                                                                                                                                                                                                     |
| 存储溶液 | PBS, 50% Glycerol                                                                                                                                                                                                          |
| 保存剂  | 0.02% Sodium azide                                                                                                                                                                                                         |
| 存放说明 | Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. |
| 浓度   | Batch dependent (Please refer to the vial label for the specific concentration.)                                                                                                                                           |
| 偶联   | Unconjugated                                                                                                                                                                                                               |
| 注意事项 | 仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。                                                                                                                                                                                 |

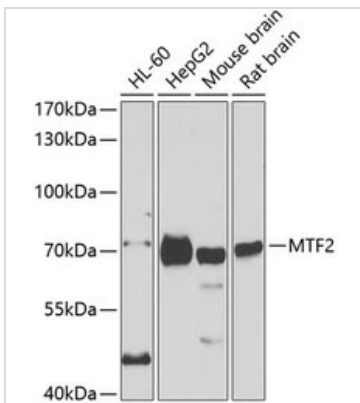


For full product information, images and publications, please visit our [website](#).

**產品圖片**

**GTX33330 ICC/IF Image**

ICC/IF analysis of A549 cells using GTX33330 MTF2 antibody.

Blue : DAPI


**GTX33330 WB Image**

WB analysis of various sample lysates using GTX33330 MTF2 antibody.

Dilution : 1:1000

Loading : 25μg per lane



For full product information, images and publications, please visit our [website](https://www.genetex.com).