

Desmin antibody [1B12]

Cat. No. GTX34350

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	WB, IHC-P
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	Assay dependent
IHC-P	1:100-1:200

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 0.5% BSA, 50% Glycerol
保存剂	0.02% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

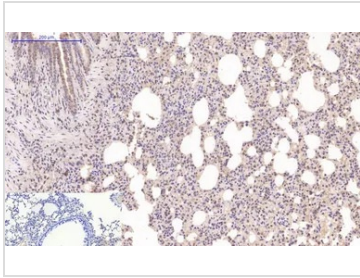


GTX34350 IHC-P Image

IHC-P analysis of mouse lung tissue using GTX34350 Desmin antibody [1B12]. Negative control (the lower left corner) was secondary antibody only.


Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

Dilution : 1:200

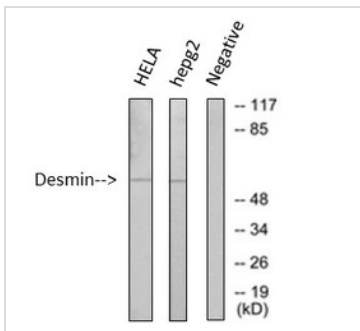

GTX34350 WB Image


WB analysis of HeLa and HepG2 cell lysates using GTX34350 Desmin antibody [1B12].

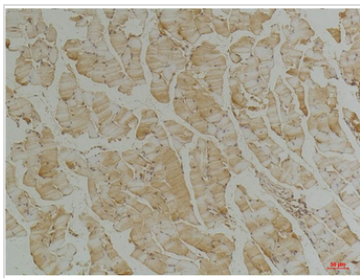

GTX34350 IHC-P Image

IHC-P analysis of mouse heart tissue using GTX34350 Desmin antibody [1B12].

Dilution : 1:200

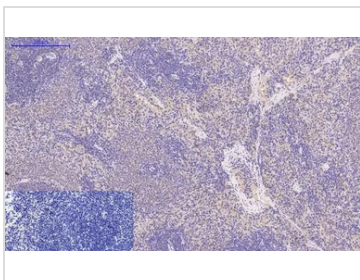


GTX34350 IHC-P Image

IHC-P analysis of mouse spleen tissue using GTX34350 Desmin antibody [1B12]. Negative control (the lower left corner) was secondary antibody only.

Antigen retrieval : Sodium citrate pH6.0 was used for antibody retrieval (>98°C, 20min)

Dilution : 1:200

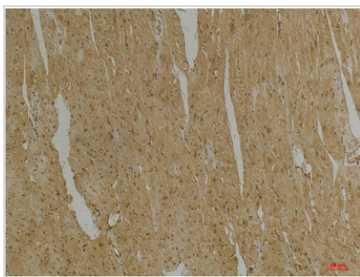

GTX34350 IHC-P Image

IHC-P analysis of mouse skeletal muscle tissue using GTX34350 Desmin antibody [1B12].

Dilution : 1:200



For full product information, images and publications, please visit our [website](http://www.genetex.com).