

SOX2 antibody [SOX2/1792]

Cat. No. GTX35087

宿主	Mouse
克隆	Monoclonal
同种型	IgG2b
实验应用	WB, IHC-P, ELISA, Protein Array
种属反应	Human, Mouse

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1-2µg/ml
IHC-P	1-2µg/ml for 30 min at RT
ELISA	2-4µg/ml (for coating)
Protein Array	Assay dependent

Note : Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes.

For ELISA coating, recommend using BSA-free format (please contact us for PBS only format).

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 0.05% BSA
保存剂	0.05% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。




For full product information, images and publications, please visit our [website](#).

產品圖片

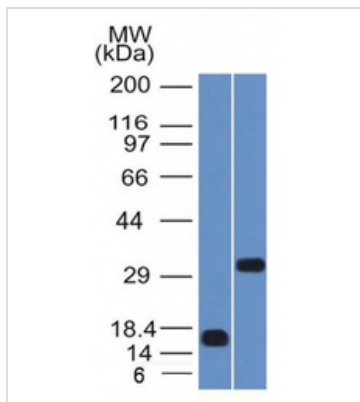


GTX35087 WB Image

WB analysis of recombinant SOX2 protein fragment and NCCIT cell lysates using GTX35087 SOX2 antibody [SOX2/1792].

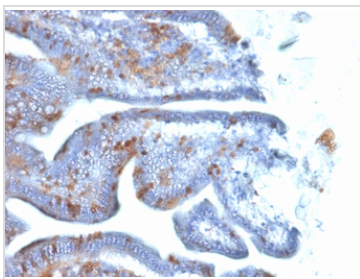


GTX35087 IHC-P Image

IHC-P analysis of mouse intestine tissue using GTX35087 SOX2 antibody [SOX2/1792].

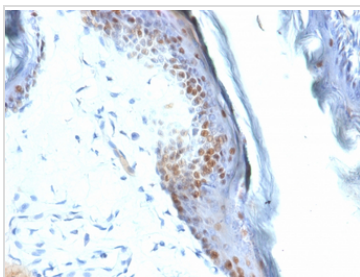

GTX35087 IHC-P Image

IHC-P analysis of mouse stomach tissue using GTX35087 SOX2 antibody [SOX2/1792].

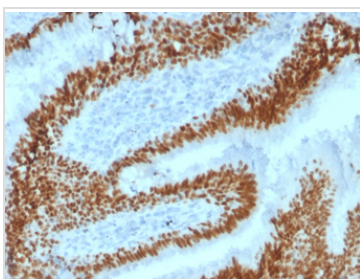

GTX35087 IHC-P Image

IHC-P analysis of human lung carcinoma tissue using GTX35087 SOX2 antibody [SOX2/1792].



For full product information, images and publications, please visit our [website](https://www.genetex.com).