

Collagen X antibody

Cat. No. GTX37732

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P, IHC
种属反应	Human, Mouse, Rat, Rabbit

引用文献 (14)

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:200-2000
IHC-P	1:300-500
IHC	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	PBS
保存剂	0.02% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.5 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_11171369
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。


For full product information, images and publications, please visit our [website](#).

產品圖片

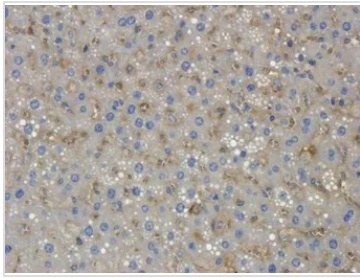


GTx37732 IHC-P Image

IHC-P analysis of rat liver tissue using GTx37732 Collagen X antibody.
Dilution : 2.5µg/ml

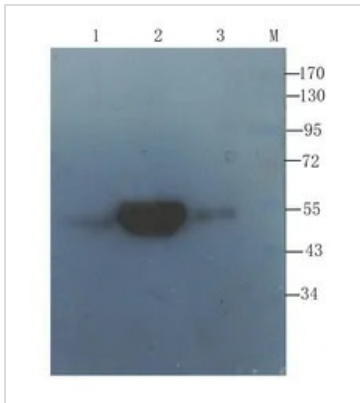


GTx37732 WB Image

WB analysis of mouse pancreas (Lane 1), human breast cancer (Lane 2), human mammary fibroma (Lane 3) lysates using GTx37732 Collagen X antibody.
Dilution : 1µg/ml

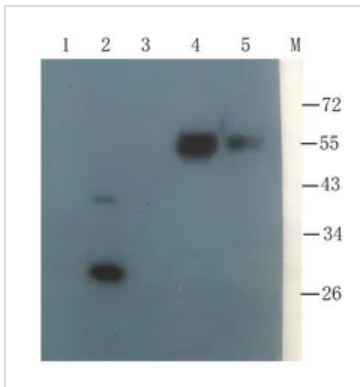

GTx37732 WB Image

WB analysis of mouse spinal cord (Lane 1), mouse liver (Lane 2), mouse pancreas (Lane 3), human breast cancer (Lane 4), human mammary fibroma (Lane 5) lysates using GTx37732 Collagen X antibody.
Dilution : 1µg/ml

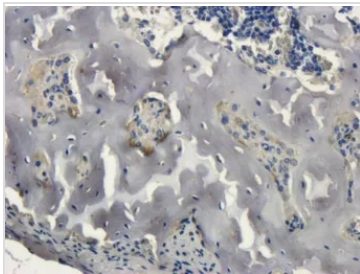

GTx37732 IHC-P Image

IHC-P analysis of mouse bone marrow tissue using GTx37732 Collagen X antibody.
Dilution : 2.5µg/ml



For full product information, images and publications, please visit our [website](https://www.genetex.com).