

TRPV5 antibody

Cat. No. GTX54868

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, ICC/IF, IHC-P
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	Assay dependent
ICC/IF	Assay dependent
IHC-P	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 1% BSA
保存剂	0.05% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.8 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).


產品圖片

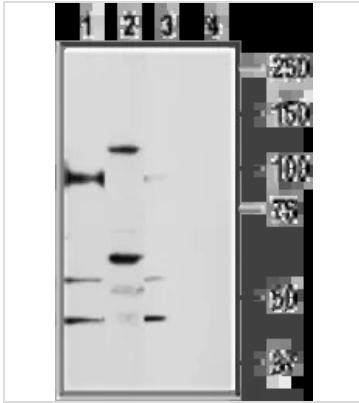

GTX54868 WB Image

WB analysis of rat kidney (lanes 1 and 3) and PC3 cell line (lanes 2 and 4) lysates using GTX54868 TRPV5 antibody preincubated with or without immunogen peptide.

Dilution : 1:200

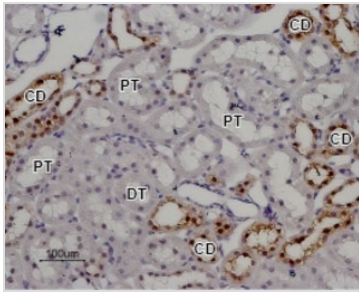


GTX54868 IHC-P Image

IHC-P analysis of rat kidney tissue using GTX54868 TRPV5 antibody. TRPV5 is expressed in collecting ducts (CD) in the renal cortex. Note that both proximal tubules (PT) and distal tubules (DT) are negative.

Hematoxylin is used as the counterstain.

Dilution : 1:100

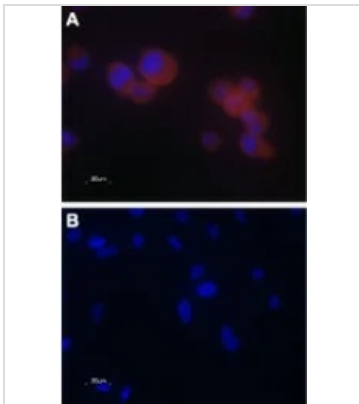

GTX54868 ICC/IF Image

ICC/IF analysis of PFA-fixed PC-3 cells using GTX54868 TRPV5 antibody.

Panel A : Primary antibody

Panel B : Primary antibody prior incubated with immunogen peptide

Dilution : 1:20



For full product information, images and publications, please visit our [website](https://www.genetex.com).