

## AHA-1 antibody [AT3E9]

Cat. No. GTX57644

|      |            |
|------|------------|
| 宿主   | Mouse      |
| 克隆   | Monoclonal |
| 同种型  | IgG2b      |
| 实验应用 | WB, ICC/IF |
| 种属反应 | Human      |

## 实验应用

## 应用说明

\*最佳稀释倍数与浓度应由研究人员确认

| Suggested dilution | Recommended dilution                     |
|--------------------|------------------------------------------|
| WB                 | Recommended starting dilution is 1:1000. |
| ICC/IF             | Assay dependent                          |

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


## 属性

|      |                                                                                                                                                                                                                            |
|------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 形式   | Liquid                                                                                                                                                                                                                     |
| 存储溶液 | PBS, 10% Glycerol                                                                                                                                                                                                          |
| 保存剂  | 0.02% Sodium azide                                                                                                                                                                                                         |
| 存放说明 | Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. |
| 浓度   | 1 mg/ml (Please refer to the vial label for the specific concentration.)                                                                                                                                                   |
| 偶联   | Unconjugated                                                                                                                                                                                                               |
| 注意事项 | 仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。                                                                                                                                                                                 |



For full product information, images and publications, please visit our [website](#).

## 產品圖片

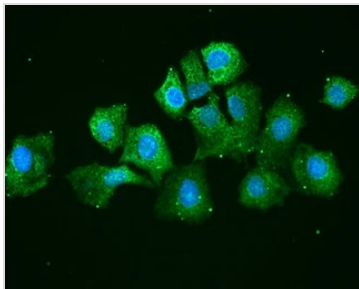


### GTX57644 ICC/IF Image

ICC/IF analysis of Hep3B cells using GTX57644 AHA-1 antibody.

Blue: DAPI

Green: Primary antibody

Dilution: 1:100

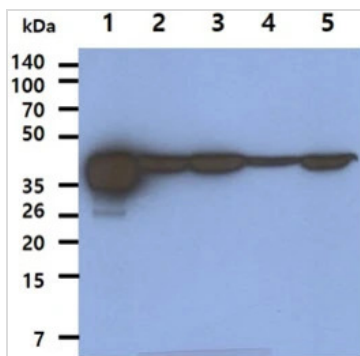

### GTX57644 WB Image

WB analysis of various samples using GTX57644 AHA-1 antibody.


Lane 1 : Recombinant Human AHA1 protein (50 ng)

Lane 2 : MCF-7 whole cell lysate (40 µg)

Lane 3 : PC-3 whole cell lysate (40 µg)

Lane 4 : 293T whole cell lysate (40 µg)

Lane 5 : HeLa whole cell lysate (40 µg)

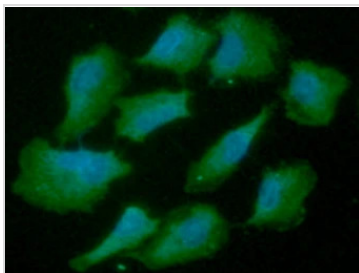

### GTX57644 ICC/IF Image

ICC/IF analysis of A549 cells using GTX57644 AHA-1 antibody.

Blue: DAPI

Green: Primary antibody

Dilution: 1:100



For full product information, images and publications, please visit our [website](#).