

GluR4 antibody

Cat. No. GTX57743

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-Fr
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	Assay dependent
IHC-Fr	Assay dependent

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 1% BSA
保存剂	0.05% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.6 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

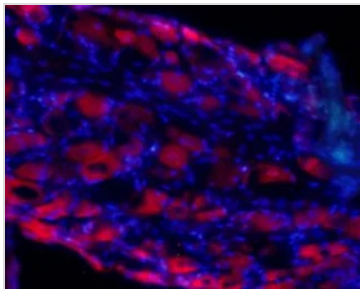


GTX57743 IHC-Fr Image

IHC-Fr analysis of rat dorsal root ganglion (DRG) tissue using GTX57743 GluR4 antibody.

Red : Primary antibody

Blue : DAPI

Dilution : 1:100

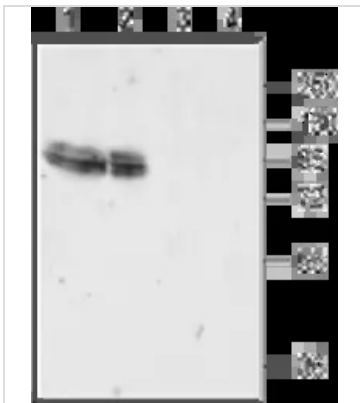

GTX57743 WB Image

WB analysis of rat brain membranes (Lanes 1 and 3) and mouse brain (Lanes 2 and 4) lysates using GTX57743 GluR4 antibody.

Lane 1 and 2 : Primary antibody

Lane 3 and 4 : Primary antibody preincubated with the control peptide antigen

Dilution : 1:200



For full product information, images and publications, please visit our [website](#).