

SOX2 antibody [10F10]

Cat. No. GTX60381

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	WB, ICC/IF, IHC-P, ELISA
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1/500 - 1/2000
ICC/IF	1/200 - 1/1000
IHC-P	1/200 - 1/1000
ELISA	1/10000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	Ascites
保存剂	0.03% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

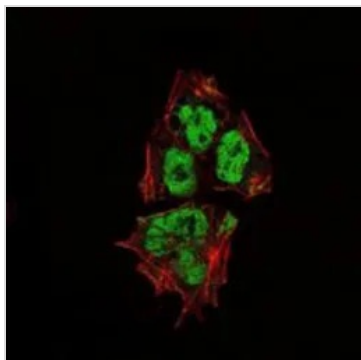


GTX60381 ICC/IF Image

ICC/IF analysis of NTERA-2 cells using GTX60381 SOX2 antibody [10F10].


Green : SOX2

Red: Actin filaments

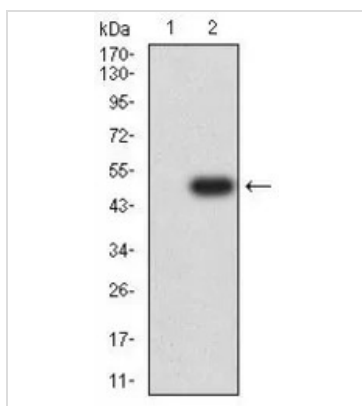

GTX60381 WB Image

WB analysis of HEK293 (1) and SOX2(AA: 2-317)-hlgGFc transfected HEK293 (2) cell lysate using GTX60381 SOX2 antibody [10F10].

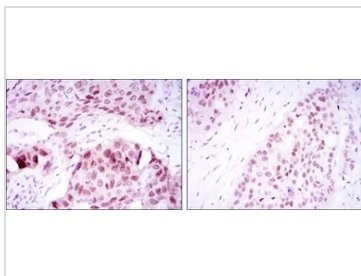

GTX60381 IHC-P Image

IHC-P analysis of lung cancer tissue (left) and esophageal cancer tissue (right) using GTX60381 SOX2 antibody [10F10].



For full product information, images and publications, please visit our [website](https://www.genetex.com).