

CRK antibody [3G11E8]

Cat. No. GTX60488

宿主	Mouse
克隆	Monoclonal
同种型	IgG2b
实验应用	WB, ICC/IF, IHC-P, FCM, ELISA
种属反应	Human, Mouse

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1/500 - 1/2000
ICC/IF	1/200 - 1/1000
IHC-P	1/200 - 1/1000
FCM	1/200 - 1/400
ELISA	1/10000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	Ascites
保存剂	0.03% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。


For full product information, images and publications, please visit our [website](#).

產品圖片

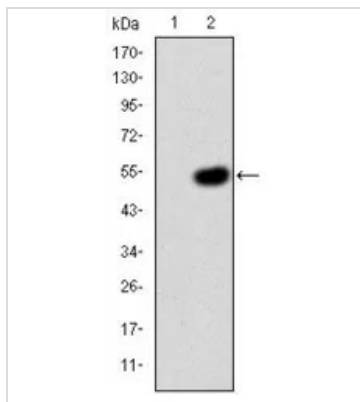

GTx60488 WB Image

WB analysis of HEK293 (1) and CRK(AA: 1-204)-hlgGfC transfected HEK293 (2) cell lysate using GTx60488 CRK antibody [3G11E8].

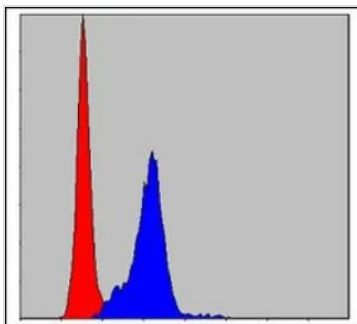


GTx60488 FCM Image

FACS analysis of HeLa cells using GTx60488 CRK antibody [3G11E8].

Blue : CRK

Red : negative control

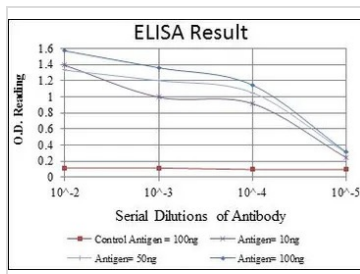

GTx60488 ELISA Image


ELISA analysis of antigen using GTx60488 CRK antibody [3G11E8].

Red : Control antigen 100ng

Purple : Antigen 10ng

Green : Antigen 50ng

Blue : Antigen 100ng

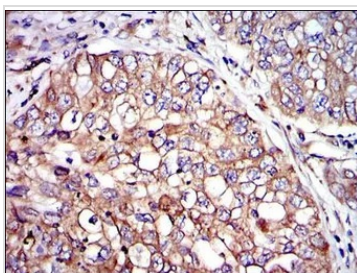

GTx60488 IHC-P Image

IHC-P analysis of lung cancer tissue using GTx60488 CRK antibody [3G11E8].



For full product information, images and publications, please visit our [website](http://www.genetex.com).