

Rabbit Anti-Human IgA antibody [RM128]

Cat. No. GTX60873

宿主	Rabbit
克隆	Monoclonal
同种型	IgG
实验应用	WB, ICC/IF, ELISA, IHC, Sandwich ELISA
种属反应	Human

Package
100 µg

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	Assay dependent
ICC/IF	0.5µg/mL - 2µg/mL
ELISA	50ng/well - 200ng/well (for Capture) ; 0.05µg/mL - 0.2µg/mL (for Detection)
IHC	0.1µg/mL - 1µg/mL
Sandwich ELISA	Assay dependent

Note : Capture : GTX60873, Detection : GTX60920

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

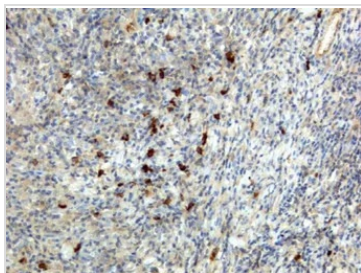
IHC-Fr: 冰冻切片

产品说明	This antibody reacts to human IgA, including both IgA1 and IgA2. No cross reactivity with human IgG, IgM, IgD, or IgE.
------	--

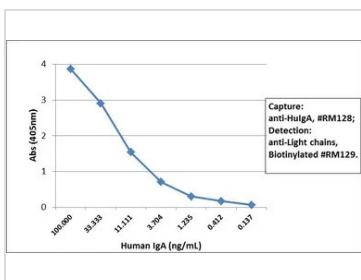
属性	
形式	Liquid
存储溶液	PBS, 1% BSA, 50% Glycerol
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
免疫原	Human IgA

For full product information, images and publications, please visit our [website](#).

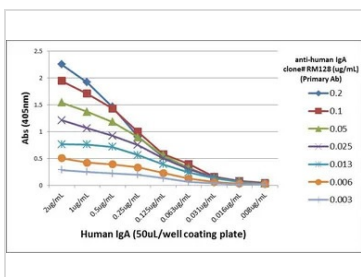
纯化方式	Protein A purified From tissue culture supernatant
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

產品图片

GTX60873 IHC Image

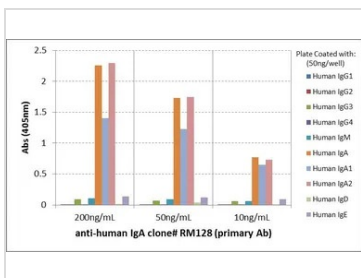
IHC analysis of human lymphoid tissue using GTX60873 Human IgA antibody [RM128].


GTX60873 ELISA Image

Sandwich ELISA analysis of human IgA using GTX60873 Human IgA antibody [RM128] as the capture antibody (100ng/well) and GTX60920 Human Ig Light Chains antibody [RM129] (Biotin) as the detection antibody.


GTX60873 ELISA Image

ELISA analysis of different amount of immunogen proteins using serial diluted GTX60873 Human IgA antibody [RM128].


GTX60873 ELISA Image

ELISA analysis of human immunoglobulins (50 ng/well) using GTX60873 Human IgA antibody [RM128]. No cross reactivity with Human IgG, IgM, IgD, or IgE.

Dilution : 200 ng/mL, 50 ng/mL, or 10 ng/mL of RM128 was used as the primary antibody.



For full product information, images and publications, please visit our [website](http://www.genetex.com).