

Opsin Red / Green antibody [HL3484]

Cat. No. GTX641369

宿主	Rabbit
克隆	Monoclonal
同种型	IgG
实验应用	IHC-P
种属反应	Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IHC-P	1:100-1:1000

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS
保存剂	No preservative
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

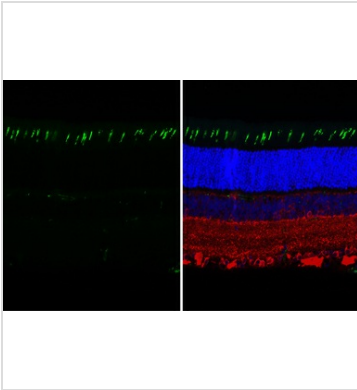

GTX641369 IHC-P Image

Green sensitive opsin antibody [HL3484] detects Green sensitive opsin protein by immunohistochemical analysis.


Sample: Paraffin-embedded rat eye.

Green: Green sensitive opsin stained by Green sensitive opsin antibody [HL3484] (GTX641369) diluted at 1:250.

Red: beta Tubulin 3/ Tuj1 stained by beta Tubulin 3/ Tuj1 antibody [GT11710] (GTX631836) diluted at 1:500.

Blue: Fluoroshield with DAPI (GTX30920).

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min

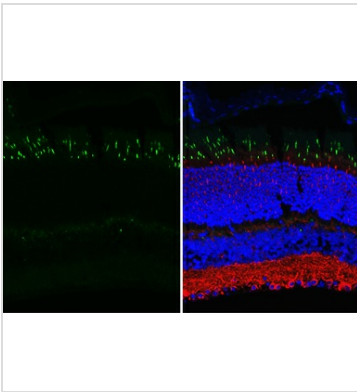

GTX641369 IHC-P Image

Green sensitive opsin antibody [HL3484] detects Green sensitive opsin protein by immunohistochemical analysis.

Sample: Paraffin-embedded mouse eye.

Green: Green sensitive opsin stained by Green sensitive opsin antibody [HL3484] (GTX641369) diluted at 1:250.

Red: beta Tubulin 3/ Tuj1 stained by beta Tubulin 3/ Tuj1 antibody [GT11710] (GTX631836) diluted at 1:500.

Blue: Fluoroshield with DAPI (GTX30920).

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min



For full product information, images and publications, please visit our [website](https://www.genetex.com).