

## ApoER2 antibody

Cat. No. GTX66733

|      |                   |
|------|-------------------|
| 宿主   | Rabbit            |
| 克隆   | Polyclonal        |
| 同种型  | IgG               |
| 实验应用 | WB                |
| 种属反应 | Human, Mouse, Rat |

## 实验应用

## 应用说明

\*最佳稀释倍数与浓度应由研究人员确认

| Suggested dilution | Recommended dilution |
|--------------------|----------------------|
| WB                 | 1:500-1:1000         |

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


## 属性

|      |                                                                                                                                                                                                                            |
|------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 形式   | Liquid                                                                                                                                                                                                                     |
| 存储溶液 | PBS                                                                                                                                                                                                                        |
| 保存剂  | 0.05% Sodium azide                                                                                                                                                                                                         |
| 存放说明 | Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. |
| 浓度   | 1 mg/ml (Please refer to the vial label for the specific concentration.)                                                                                                                                                   |
| 纯度   | > 95% (by SDS-PAGE)                                                                                                                                                                                                        |
| 偶联   | Unconjugated                                                                                                                                                                                                               |
| 注意事项 | For In vitro laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.                                                                          |



For full product information, images and publications, please visit our [website](#).

產品圖片

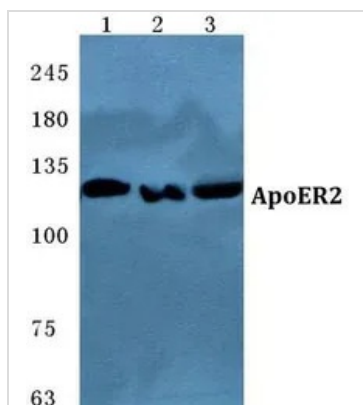

**GTX66733 WB Image**

WB analysis of various samples using GTX66733 ApoER2 antibody.

Lane1 : MCF-7 cell lysate

Lane2 : NIH-3T3 cell lysate

Lane3 : H9C2 cell lysate

Dilution : 1:500



For full product information, images and publications, please visit our [website](https://www.genetex.com).