

## mGluR6 antibody

Cat. No. GTX71531

|      |                                                                   |
|------|-------------------------------------------------------------------|
| 宿主   | Rabbit                                                            |
| 克隆   | Polyclonal                                                        |
| 同种型  | IgG                                                               |
| 实验应用 | WB, IHC-P                                                         |
| 种属反应 | Human, Mouse, Rat, Rabbit, Bovine, Dog, Hamster, Pig, Monkey, Bat |

## 实验应用

## 应用说明

\*最佳稀释倍数与浓度应由研究人员确认

| Suggested dilution | Recommended dilution |
|--------------------|----------------------|
| WB                 | Assay dependent      |
| IHC-P              | 12 µg/ml             |

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


## 属性

|      |                                                                                                                                                                                                                            |
|------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 形式   | Liquid                                                                                                                                                                                                                     |
| 存储溶液 | PBS                                                                                                                                                                                                                        |
| 保存剂  | 0.1% Sodium azide                                                                                                                                                                                                          |
| 存放说明 | Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. |
| 浓度   | 1 mg/ml (Please refer to the vial label for the specific concentration.)                                                                                                                                                   |
| 偶联   | Unconjugated                                                                                                                                                                                                               |
| RRID | AB_375799                                                                                                                                                                                                                  |
| 注意事项 | 仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。                                                                                                                                                                                |




For full product information, images and publications, please visit our [website](#).

## 產品圖片

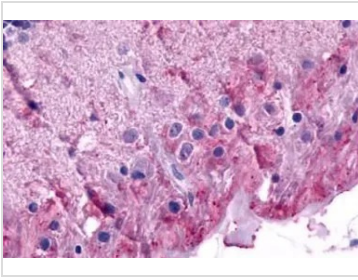

### GTX71531 IHC-P Image


IHC-P analysis of eye, retina, ganglion cells tissue using GTX71531 mGluR6 antibody.

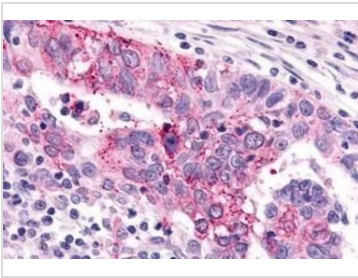

### GTX71531 IHC-P Image

IHC-P analysis of human ovary, carcinoma tissue using GTX71531 mGluR6 antibody.


Antigen retrieval : Heat-induced antigen retrieval

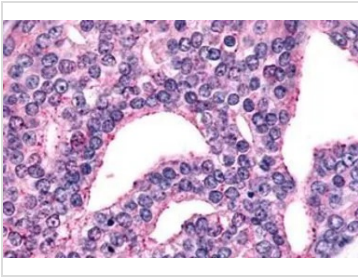

### GTX71531 IHC-P Image

IHC-P analysis of human prostate, carcinoma tissue using GTX71531 mGluR6 antibody.


Antigen retrieval : Heat-induced antigen retrieval

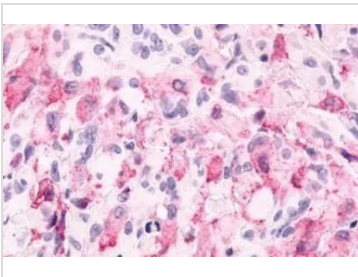

### GTX71531 IHC-P Image

IHC-P analysis of human brain, glioblastoma tissue using GTX71531 mGluR6 antibody.

Antigen retrieval : Heat-induced antigen retrieval

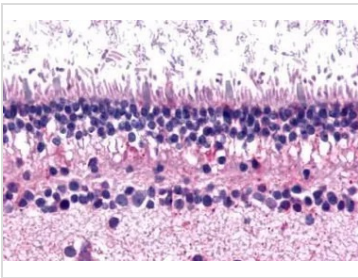

### GTX71531 IHC-P Image

IHC-P analysis of eye, retina tissue using GTX71531 mGluR6 antibody.



For full product information, images and publications, please visit our [website](#).