

Integrin alpha 4 antibody [44H6]

Cat. No. GTX76621

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	IHC-Fr, FCM
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
IHC-Fr	1/50-1/100
FCM	1/40

Note : Overnight incubation is recommended. The epitope recognised by this antibody is reported to be sensitive to formaldehyde fixation and tissue processing. We recommends the use of acetone fixation for frozen sections.

Use 10µl of the suggested working dilution to label 100µl of human whole peripheral blood.

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

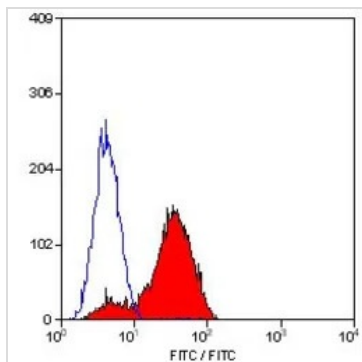
IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	PBS
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	1.0 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
RRID	AB_377460
注意事项	仅供实验室使用。不适用于人类或动物的任何临床, 治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

GTX76621 FCM Image

FACS analysis of human peripheral blood lymphocytes using GTX76621 Integrin alpha 4 antibody [44H6].



For full product information, images and publications, please visit our [website](#).