

Goat Anti-Dog IgG antibody (HRP)

Cat. No. GTX77014

宿主	Goat
克隆	Polyclonal
同种型	IgG
实验应用	WB, ELISA, IHC
种属反应	Dog

Package
500 µg

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:1000-1:10000
ELISA	1:50000-1:150000
IHC	1:500-1:2500

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	20mM Potassium Phosphate, 150mM NaCl, 1% BSA
保存剂	0.01% Gentamicin Sulfate
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	2 mg/ml (Please refer to the vial label for the specific concentration.)
免疫原	Dog IgG whole molecule
纯化方式	Purified by antigen-affinity chromatography. From serum
偶联	Horseradish peroxidase(HRP)
注意事项	仅供实验室使用。不适用于人类或动物的任何临床，治疗或诊断用途。不适合动物或人类食用。

For full product information, images and publications, please visit our [website](#).

產品圖片

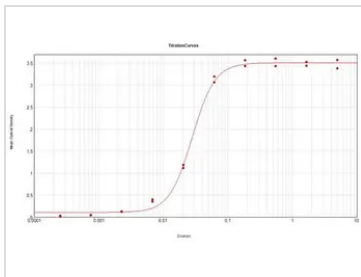

GTx77014 ELISA Image

ELISA analysis of dog IgG using serially diluted GTx77014 Goat Anti-Dog IgG antibody (HRP).
Coating : 1 µg



For full product information, images and publications, please visit our [website](https://www.genetex.com).