

COMP antibody, Internal

Cat. No. GTX81660

宿主	Rabbit
克隆	Polyclonal
同种型	IgG
实验应用	WB, IHC-P, FCM
种属反应	Human, Mouse, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:2000
IHC-P	1:10-1:50
FCM	1:25

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	Batch dependent (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

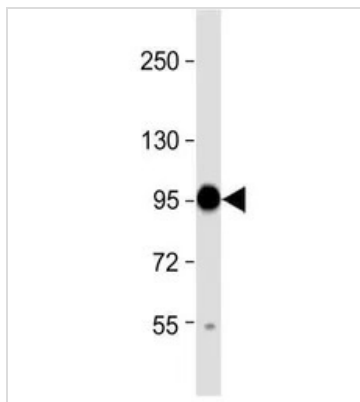


GTx81660 WB Image

WB analysis of Rat skeletal muscle lysate using GTX81660 COMP antibody, Internal.

Loading : 20 µg per lane

Dilution : 1:2000

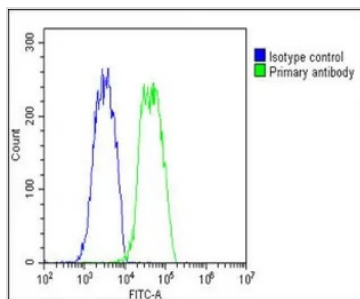

GTx81660 FCM Image

FACS analysis of HepG2 cells using GTX81660 COMP antibody, Internal.


Green : primary antibody

Blue : Isotype control

Fixation : 2% paraformaldehyde (PFA) (10 min)

Permeabilization : 90% methanol (10 min)

Dilution : 1:25

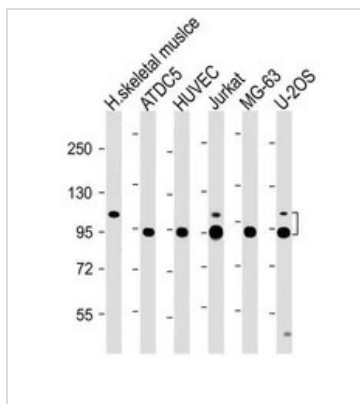

GTx81660 WB Image

WB analysis of various samples using GTX81660 COMP antibody, Internal.

Lane 1: Human skeletal muscle lysate

Lane 2: ATDC5 whole cell lysate

Lane 3: HUVEC whole cell lysate


Lane 4: Jurkat whole cell lysate

Lane 5: MG-63 whole cell lysate

Lane 6: U-2OS whole cell lysate

Loading : 20 µg per lane

Dilution : 1:2000

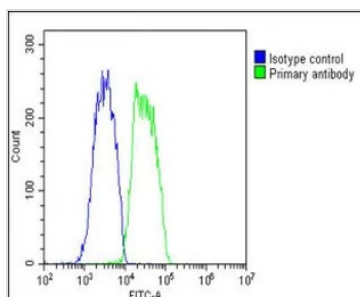

GTx81660 FCM Image

FACS analysis of HepG2 cells using GTX81660 COMP antibody, Internal.

Green : primary antibody

Blue : Isotype control

Fixation : 2% paraformaldehyde (PFA) (10 min)

Permeabilization : 90% methanol (10 min)

Dilution : 1:25



For full product information, images and publications, please visit our [website](https://www.genetex.com).