

PPME1 antibody [7F12]

Cat. No. GTX83847

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	WB, ICC/IF, IHC-P, FCM
种属反应	Human, Mouse, Rat, Dog, Monkey

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	1:2000
ICC/IF	1:100
IHC-P	1:50
FCM	1:100

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 1% BSA, 50% Glycerol
保存剂	0.02% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.96 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

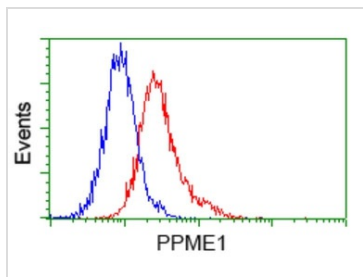


GTX83847 FCM Image

FACS analysis of Jurkat cells using GTX83847 PPME1 antibody [7F12].

Red : Primary antibody

Blue : Negative control antibody

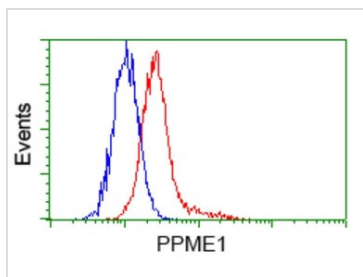


GTX83847 FCM Image

FACS analysis of HeLa cells using GTX83847 PPME1 antibody [7F12].


Red : Primary antibody

Blue : Negative control antibody

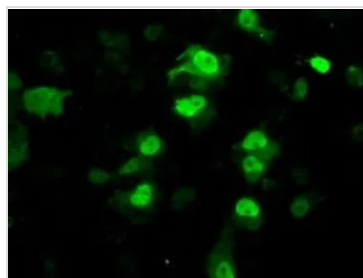

GTX83847 ICC/IF Image

ICC/IF analysis of COS7 cells transiently transfected with PPME1 plasmid using GTX83847 PPME1 antibody [7F12].

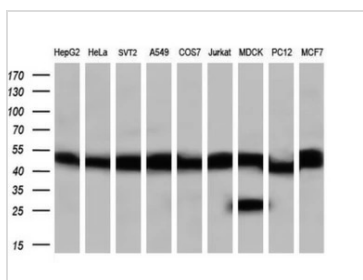

GTX83847 WB Image

WB analysis of various cell lines using GTX83847 PPME1 antibody [7F12].

Loading : 35 ug per lane

Dilution : 1:200



For full product information, images and publications, please visit our [website](#).