

CMG1 antibody, C-term

Cat. No. GTX89771

宿主	Goat
克隆	Polyclonal
同种型	IgG
实验应用	WB, ICC/IF, IHC-P
种属反应	Human, Rat

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
WB	0.3-1µg/ml
ICC/IF	Assay dependent
IHC-P	5-10µg/ml

Note : Human Liver shows vesicular staining in cytoplasm of hepatocytes.

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片


IHC-Fr: 冰冻切片

属性

形式	Liquid
存储溶液	TBS, 0.5% BSA
保存剂	0.02% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
浓度	0.50 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Unconjugated
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。

For full product information, images and publications, please visit our [website](#).

產品圖片

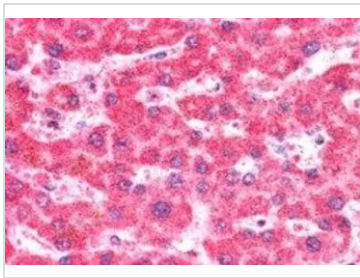


GTX89771 IHC-P Image

IHC-P analysis of human liver using GTX89771 CMG1 antibody, C-term.

Antigen retrieval : citrate buffer pH 6

Dilution : 5μg/ml

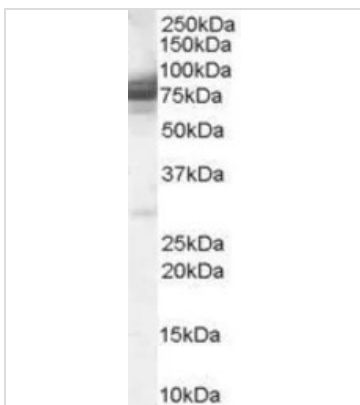


GTX89771 WB Image

WB analysis of 293 lysate using GTX89771 CMG1 antibody, C-term.

Dilution : 0.3μg/ml

Loading : 35μg protein in RIPA buffer

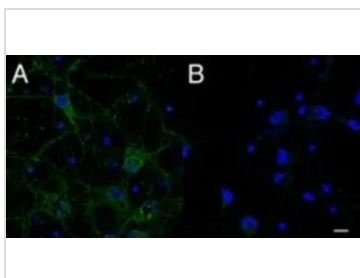

GTX89771 ICC/IF Image

ICC/IF analysis of primary rat cortical neurons using GTX89771 CMG1 antibody, C-term.

Dilution : 2.5μg/ml

A: Primary antibody, B: control



For full product information, images and publications, please visit our [website](#).